


Experiment No. 4:-

Determine the amount of Bi^{+3} gravimetrically.

Chemical Equation:-



Precipitating Agent:-

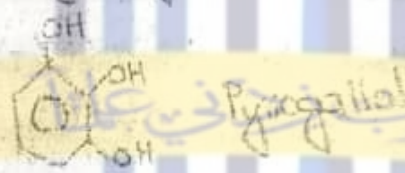

Join our what's app groups
Umer Abdullah GPGC
contact no: 0334-2194600

Observations & calculations:-

wt. of filter paper = w_1 = 0.9g

wt. of filter paper + precipitates =

w_2 = 1.2g

wt. of precipitates = w_3 = $1.2 - 0.9$

= 0.3g

3.3g of complex contain Bi^{+3} = 2.8g

1g of complex contain Bi^{+3} = $\frac{2.8}{3.3}$

0.3g of complex contain Bi^{+3} = $\frac{2.8}{3.3} \times 0.3$

= 0.25g

Experiment No. 4:-


Determine the amount/L of Bi^{+3} gravimetrically.

Principle:-

it is a gravimetric analysis.

Chemistry of Bi^{+3} :-

Bismuth is the heaviest member of all naturally occurring elements. its atomic no. is 83 and atomic wt. is 208 g/mol. It is post transition metal. it belongs to the 6th period of the periodic table. the electronic configuration of bismuth is given as follows in the below:-

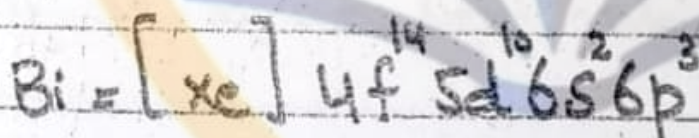

Bismuth is used in some pharmaceuticals. The medicine of Bismuth is Bismol which is used for acidity. Bismuth shows s-inert pair effect. it shows variable oxidation states like as Bi^{+3} & Bi^{+5} , it is toxic.

Procedure:-

First of all 1% of $\text{Bi}(\text{NO}_3)_3 \cdot 5\text{H}_2\text{O}$

11

Amount of Sample contain $B_{i+3} = 0.1g$

1ml of Sample contain $B_{i+3} = \frac{0.1}{20}$

Amount of Sample contain $B_{i+3} = \frac{0.1}{20} \times 1000$

$= 5g/L$

Result:-

The amount/L of B_{i+3} is 9.4g/L

PUACP

Material provided by: ece aylin
Arranged by: Umer Abdullah (UA)

11


Solution was prepared by taking 1g of $\text{Bi}(\text{NO}_3)_3 \cdot 5\text{H}_2\text{O}$ in a flask and small quantity of HNO_3 was added in it, in order to dissolve it in water. Then it was marked upto 100ml with distilled water. 1% solution of Pyrogallol was prepared by dissolving 1g of it in hot water and mark it upto 100ml with distilled water. kept the Pyrogallol solution in cap bottle, it is reactive and explosive, then 20ml of metal solution was added in 500ml beaker and Pyrogallol was added in it (7ml) till ppt. formed. If ppt. was not formed then added a few drops of liq. NH_3 and yellow coloured ppts. were formed, coagulate the ppt. by heating, filter it on dry filter paper. Then dry it in oven and note down the weight.

12


Experiment No. 6:-

Determine the amount/L of NO_3^- ions gravimetrically.

Chemical Equation:-

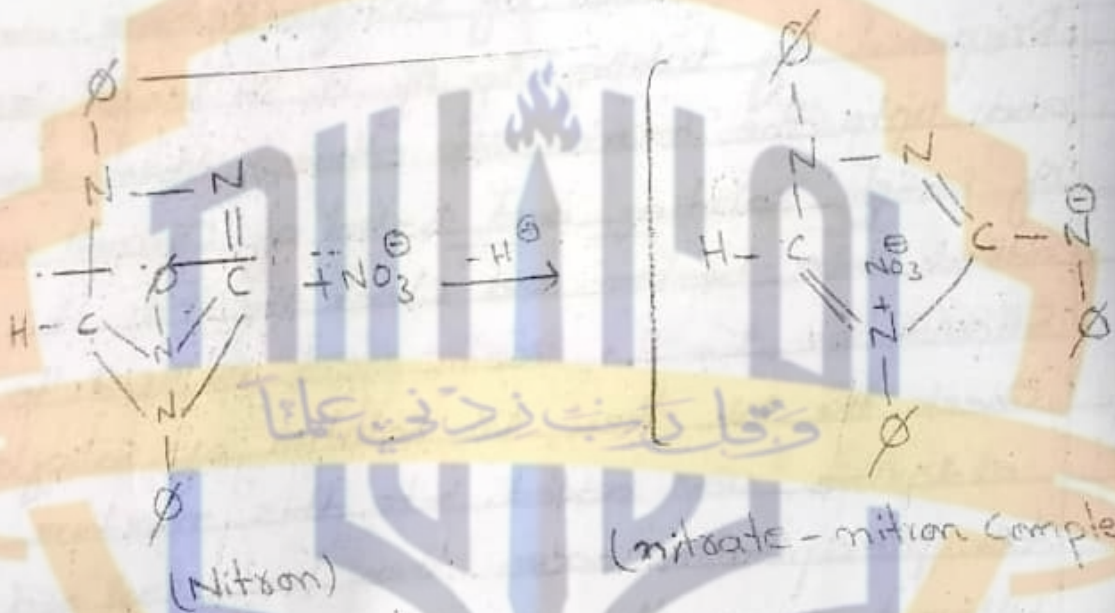

Precipitating Agent:-

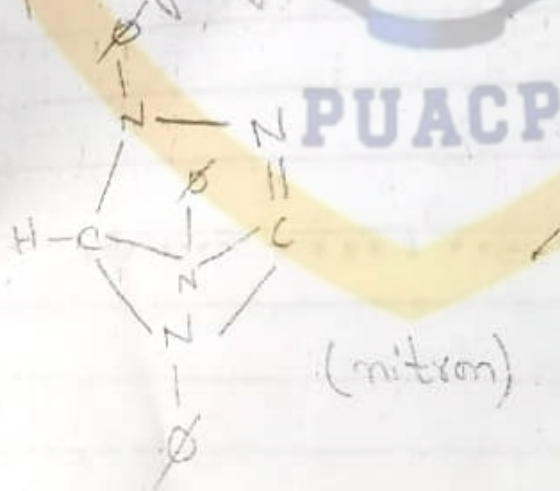

Experiment No. 6:-

Determine the amount/L of NO_3^- ion gravimetrically:-

Principle:-

It is a gravimetric analysis.

Chemistry of NO_3^- ion:-

Nitrate is polyatomic ion with molecular formula NO_3^- and a molecular weight of 62.0049 g/mol. They are mainly produced for used as fertilizers. because of their high solubility and biodegradability, the main nitrates are ammonium, potassium, sodium and calcium salts. To treat acidic soil in Pakistan fertilizers like ammonium nitrate and potassium nitrate are used. the second major application of nitrates is as an oxidizing agent.

Sodium nitrate is used to remove air bubbles from molten glass and some ceramics. they oxidizes the Fe atom in hemoglobin from ferrous Iron (+2) to ferric iron (+3) make it unable to carry O_2 . Excessive NO_3^- are runoff with rain H₂O where it causes algal bloom due to which death of the Phytoplankton and the Zooplankton.

Observations and calculations:-

$$\text{wt. of filter Paper} = w_1 g = 0.9 g$$

$$\text{wt. of Paper + ppt.} = w_2 g = 2.23 g$$

$$\text{wt. of Precipitates} = w_3 g = w_2 - w_1 = 2.23 - 0.9 = 1.33 g$$

$$374 g \text{ of complex contain } NO_3^- = 62 g$$

$$1 g \text{ of complex contain } NO_3^- = \frac{62}{374} g$$

$$1.33 g \text{ of complex contain } NO_3^- = \frac{62}{374} \times 1.33 = 0.22 g$$

$$2 \text{ ml of Sample contain } NO_3^- = 0.22 g$$

$$1 \text{ ml of Sample contain } NO_3^- = \frac{0.22}{20} g$$

$$1000 \text{ ml of Sample contain } NO_3^- = \frac{0.22}{20} \times 1000 = 11 g/l$$

Result:-

The amount of nitrate ion is $11 g/l$.




Procedure:-

First of all, 5% solution of CH_3COOH was prepared. In lab 100% acetic acid is available, so according to $C_1V_1 = C_2V_2$, 5ml of CH_3COOH was taken in a 100ml flask. 10-15ml of ethanol was added into it and upto mark with distilled water. then 1% solution of nitron was prepared by adding 1g of nitron in 100ml flask and upto mark with 5% solution of acetic acid with ethanol. if not dissolve then heat it on water bath. then 20ml of sample solution was taken in beaker. 1% solution of NO_2^- ion was prepared by dissolving 1g of NaNO_2 in 100ml flask and mark upto with distilled H_2O . then 20ml sample solution + 1ml of CH_3COOH was add in beaker, then nitron solution was added into it till precipitation occur. then Coagulate the ppt. by boiling, cool, filter it on dry filter paper and note down weight.

Experiment No. 81-

Determine amount/L of Cu^{2+} and Ni^{2+} in give mix.
by gravimetrically:-

Chemical Equations:-

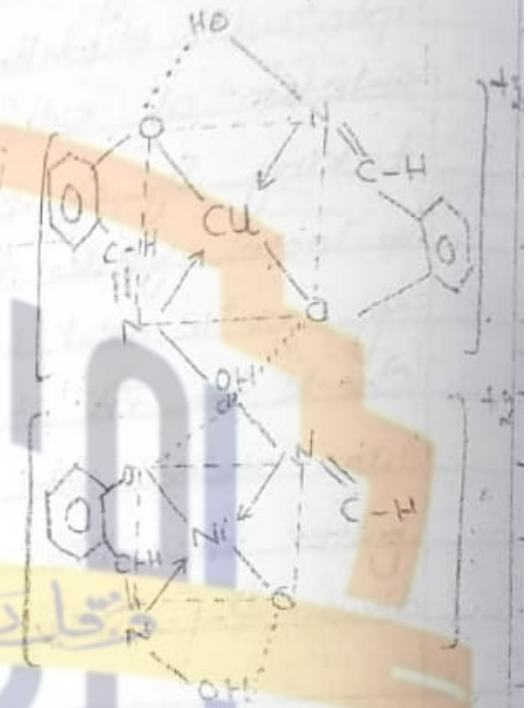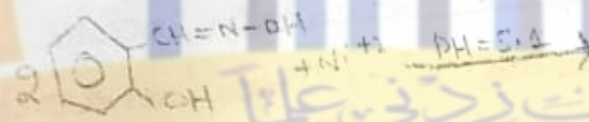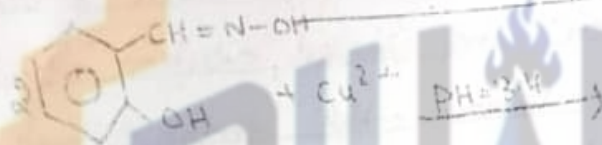

Observations and calculations:-

wt. of filter paper = w_1 g

wt. of Paper + ppt = w_2 g

wt. of Precipitate = w_3 g = $w_2 - w_1$ = Ag

PUACP

335.5g of Cu-salicylaldehyde complex contain

$$\text{Cu}^{2+} = 63.54 \text{ g}$$

1g of Cu-salicylaldehyde complex contain

$$\text{Cu}^{2+} = \frac{63.54}{335.5} \text{ g}$$

Ag of Cu-salicylaldehyde complex contain

$$\text{Cu}^{2+} = \frac{63.54}{335.5} \text{ g} \times A = 'B' \text{ g}$$

Experiment No. 8:-

Determine the amount/L of Cu^{+2} and Ni^{+2} in given mixture by gravimetrically:-

Principle:-

This is gravimetric analysis.

Theory:-

Gravimetry is non-instrumental technique and gravimetric analysis describes a set of methods in analytical chemistry for the qualitative determination of analyte based on mass.

When based on instrumental, then it is called thermal gravimetry. It is used as qualitative as well as for the quantitative analysis. Two types of errors come.

* Positive error.

* Negative error.

Positive error is due to impurity or due to deposition of organic solvents. Negative error comes when organic solvent used in large quantity and complex dissolve in it.

15

$$\begin{aligned} 20 \text{ ml of Sample contain } \text{Cu}^{2+} &= \frac{B}{20} \text{ g} \\ 1 \text{ ml of Sample contain } \text{Cu}^{2+} &= \frac{B}{20} \text{ g} \\ 1000 \text{ ml of Sample contain } \text{Cu}^{2+} &= \frac{B}{20} \times 1000 = 'C' \text{ g/L} \end{aligned}$$

$$\text{Wt. of filter paper} = W_4 \text{ g}$$

$$\text{Wt. of Paper + ppt} = W_5 \text{ g}$$

$$\text{Wt. of ppt.} = W_6 \text{ g} = W_5 - W_4 = 'D' \text{ g}$$

$$330.7 \text{ g of Ni-Salicylaldoxime complex contain } \text{Ni}^{2+} = 58.0 \text{ g}$$

$$1.00 \text{ g of Ni-Salicylaldoxime complex contain } \text{Ni}^{2+} = \frac{58}{330.7}$$

$$D \text{ g of Ni-Salicylaldoxime complex contain } \text{Ni}^{2+} = \frac{58}{330.7} \times D = 'E' \text{ g}$$

$$20 \text{ ml of Sample solution contains } \text{Ni}^{2+} = E \text{ g}$$

$$1 \text{ ml of Sample solution contains } \text{Ni}^{2+} = \frac{E}{20} \text{ g}$$

$$1000 \text{ ml of Sample solution contains } \text{Ni}^{2+} = \frac{E}{20} \times 1000 = 'F' \text{ g/L}$$

$$\text{Total amount of mixture} = G$$

$$\% \text{ age of } \text{Cu}^{2+} = \frac{C}{G} \times 100 = 'X' \%$$

$$\% \text{ age of } \text{Ni}^{2+} = \frac{F}{G} \times 100 = 'Y' \%$$

Result:-

The amount/L of Cu^{2+} is 'C' g/L and of Ni^{2+} is 'F' g/L in given mixture.

Procedure:-

0.1 g of Ni salt and 0.1 g of Cu salt was taken in 100ml flask and made upto the distilled water. it is 0.1% of sample solution of mixture. 0.1% of Salicyaldoxime in ethanol was prepared by taking the 0.1 g of it and dissolve it into 100ml of ethanol. 2ml of sample was taken and then added the fresh glacial acetic acid to maintain pH 2 or 3, then added Salicyaldoxime till the ppt. formed. for Coagulation heat the ppt, then filtered and dried it. these are Cu^{+2} ppt. then added dilute NH_3 soln. in the above filtrate and adjust pH 6 and then added the Salicyaldoxime till the precipitates formed, then heat for coagulation, filtered it and dried the ppt. in oven. these are Ni^{+2} precipitates.

Material provided by: ece ^{ay}bin
Arranged by: Umer Abdullah (UA)

Experiment no: 26

1


Determine the amount of Iron colorimetrically in given sample:-

→ Preparation of Solution:-
($\text{FeSO}_4 \cdot 7\text{H}_2\text{O}$)

$$\frac{\text{Mol. wt.}}{\text{molar mass}} = \frac{178.01}{56} = \frac{3.71}{10} = 0.371 \text{ g/100ml}$$

→ 0.25% 1,10-phenanthroline solution:-
0.25g in 100ml flask.
 $\lambda = 515\text{nm}$

→ Chemical Equations:-

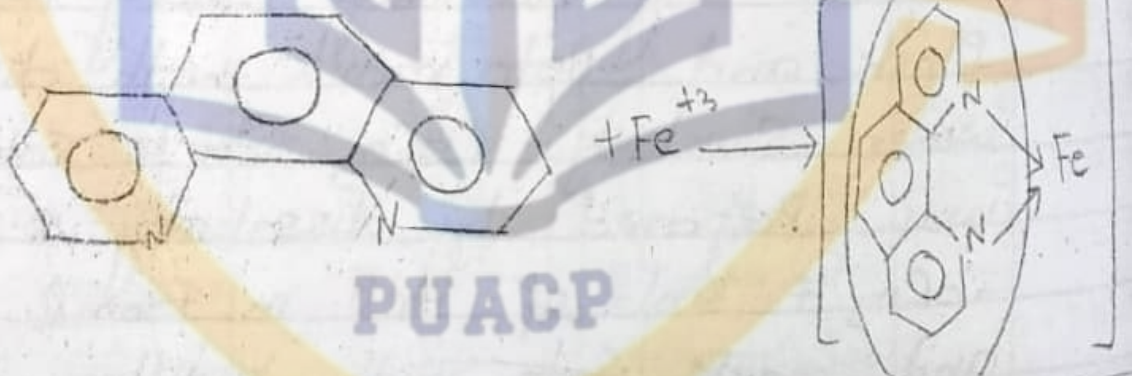

→ Buffer of PH 4.5:-

→ 0.1M CH_3COOH :-

Given : Required

$$C_1 V_1 = C_2 V_2$$

$$17 \times V_1 = 0.1 \times 100$$

$$V_1 = 0.58\text{ml/100ml}$$

→ 0.1M CH_3COONa

mol. wt = 82g

$$M = \frac{\text{mass} \times 1000}{\text{m. mass} \times 100}$$

$$0.1 = \frac{x}{100} \times 100$$

$$x = \frac{82}{0.1 \times 82}$$

$$= 0.82\text{g/100ml}$$

Experiment No : 26

Determine the amount/L of Iron Colorimetrically in given Sample:-

Principle:-

This is a colorimetric analysis of iron and fall in visible spectroscopy.

Theory:-

Spectroscopy is interaction of radiations with matter. it can be classified on the basis of wavelength.

These wavelengths are following like X-ray, U.V, visible, I.R, microwave Spectroscopy and γ -ray Spectroscopy.

On the basis of Species,
Atomic Spectroscopy.

Flame-photometry

AAS

Tcp-Analysis

Molecular Spectroscopy.

U.V

Visible

3

→ For Buffer preparation:-

Mix 65 ml of CH_3COOH and 35ml of 0.1M of CH_3COONa .

→ 7% Hydroquinone:-


1g hydroquinone in 10ml of buffer. to dissolve hydroquinone, 5-6ml of $\text{C}_2\text{H}_5\text{OH}$ was added into it.

→ For standard solution preparation:-


In a test tube 10ml of standard solution + 2ml of hydroquinone + 5ml of orthophenanthroline was added.

→ For blank solution:-

10ml of dis. H_2O + 2ml of hydroquinone + 5ml of orthophenanthroline was added into it.

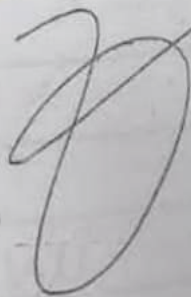

1. Visible Spectroscopy is called colorimetry.

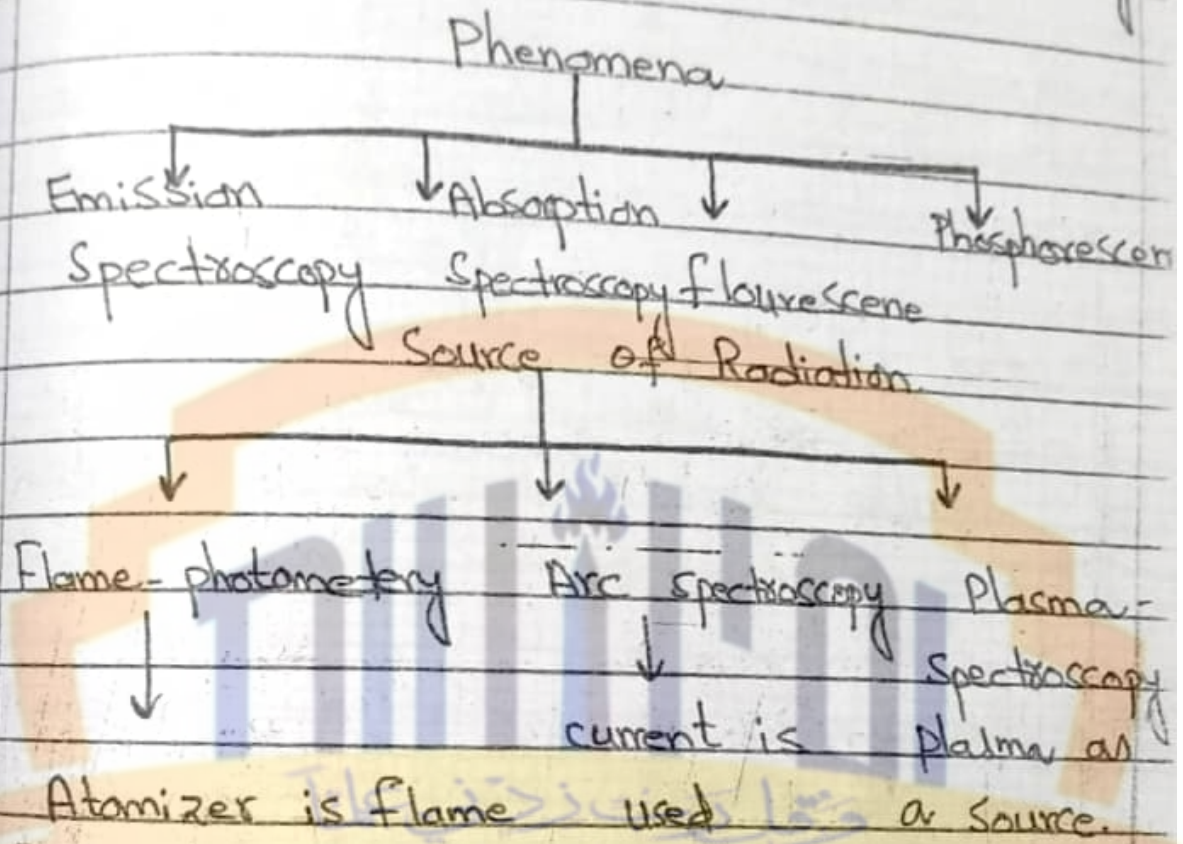

Basic Requirements of colorimetric analysis:-
 Species would be soluble in either solvent (water or organic), suspensions are not allowed.

It should be coloured e.g. in organic Caffeine or Ni-analysis give yellow colour.

Procedure:-

1% hydroquinone was prepared by dissolving 1g of hydroquinone in 100ml of buffer.
 If hydroquinone is not dissolved, then add 5-6ml of ethanol to dissolve the hydroquinone.

Scale along X

Large scale = 10


Small scale = 1.0

Scale along y

Large scale = 0.2

Small scale = 0.02

5

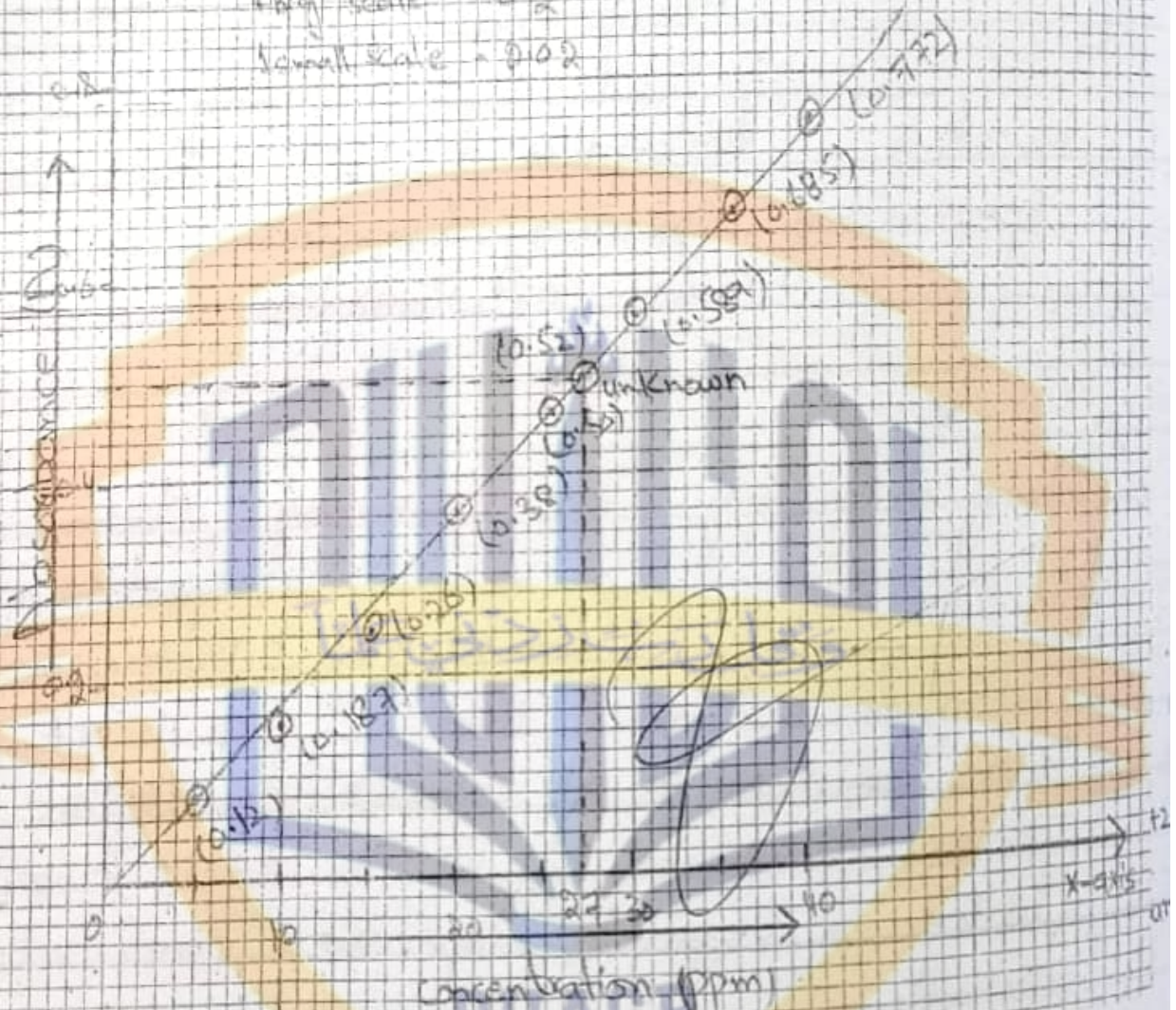

Ni 5-50ppm Solution were prepared. 5ppm solution was prepared by dissolving 1ml standard solution. 2ml of hydroquinone and 5ml of phenanthroline in test tube. Similarly other solutions were prepared.

Reading was taken with the help of the electronic. it was set at 2.0 with the help of blank solution. And blank solution was prepared by dissolving 1ml of dist. H_2O + 2ml of hydroquinone + 5ml of the orthophenanthroline. Instrument was set at 515nm and absorbance was noted, and graph was plotted b/w absorbance and concentration.

Experiment No: 27

Determine the amount/L of Ni^{2+} by colorimetrically:-

→ Preparation of Solutions:-

($\text{NiSO}_4 \cdot 6\text{H}_2\text{O}$)

$$\frac{\text{Mol. wt.}}{\text{m.mass}} = \frac{263}{58} \cdot \frac{4.53}{10} = 0.453 \text{ g/100ml}$$

→ 1% DMG solution:-

0.1g of DMG in ethanol

→ 0.5% NH_3 Soln:-

Given : Required

Material provided by: ece aydin
Arranged by: Umer Abdullah (UA)

$$C_1 V_1 = C_2 V_2$$

$$33 \times V_1 = 5 \times 100$$

$$V_1 = 15 \text{ mL in 100 mL flask}$$

→ For standard solution preparation:-

In a Separating funnel 10 mL CHCl_3 +
10 mL standard solution + 1 drop of
 NH_3 + 10 mL DMG.

Lower organic layer separate out.

Experiment No: 87

of Determine the amount/L of Ni^{+2} by
Colorimetrically in a given sample.

Principle:-

This is colorimetric analysis
of Nickel.

Theory:-

Spectroscopy is interaction of
radiations with matter. It can be
classified on the basis of wavelength
and the species. Visible Spectroscopy
is Colorimetry.

Basic requirements of Colorimetric analysis:
Species should be soluble in either
solvent. Suspensions are not to be
allowed.

it should be coloured or form coloured
complex.

Law of Colorimetric analysis is the
Beer-Lambert law;

$$A \propto CL$$


$$A = \epsilon CL$$

Procedure:-

First of all preparation of following soln.
are done.

Scale along x-axis
 1 big scale = 10
 1 small scale = 1.0

Scale along y-axis
 1 big scale = 0.2
 1 small scale = 0.02

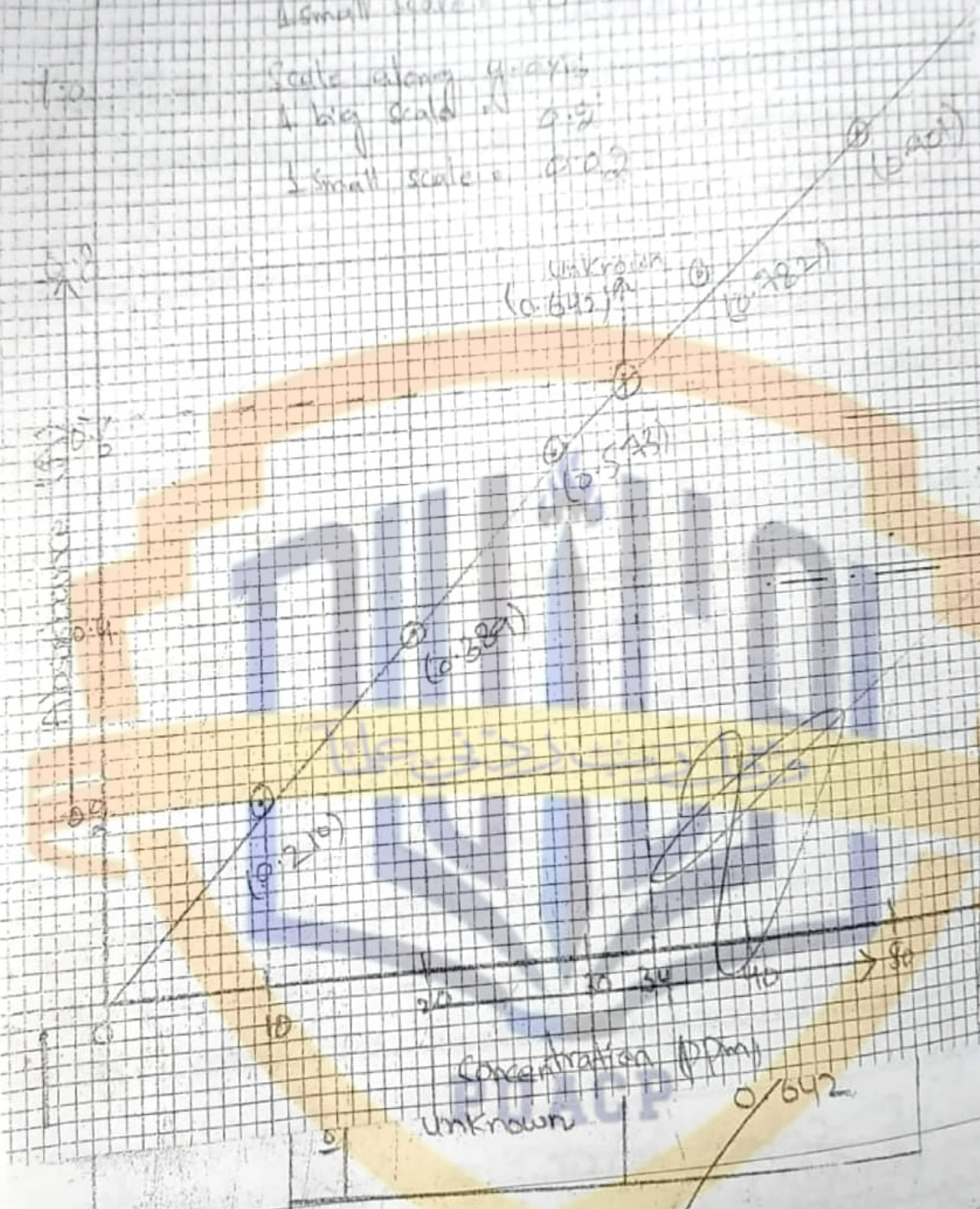

→ Result:-

The conc. of unknown sample is at absorbance of 0.642 is 34 ppm.

Metal ion solution i.e. Ni^{2+} 0.45g of $\text{NiSO}_4 \cdot 6\text{H}_2\text{O}$ in 100ml flask and upto mark it with distilled H_2O . 0.1% DMG solution was prepared in ethanol. 5% NH_3 solution was prepared according to dilution formula.

$$C_1 V_1 = C_2 V_2$$

15 ml of NH_3 dissolved in 100ml flask. In separating funnel, 10ml of CHCl_3 , 10ml of standard solution and 1 drop of NH_3 and 10ml of DMG, was added in the flask. Shaked and separate the lower organic layer, which is of yellow colour.



Spectronic was set at 0 by putting the chloroform as blank solution.

Readings were taken with the help of Spectronic 20. The graph was plotted of concentration and absorbance which is straight line.

Experiment No: 28

Determine the amount/L of Cr^{+6} by using diphenyl carbazide colorimetrically.

→ Chemical Equation:-

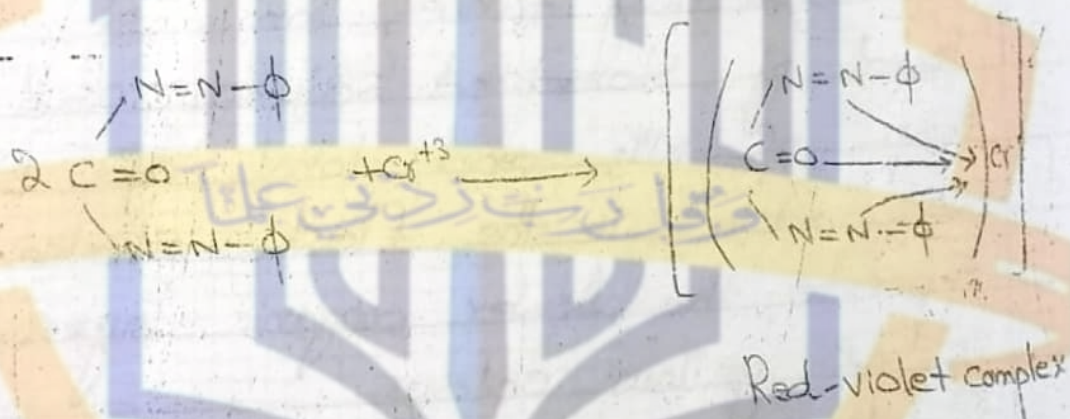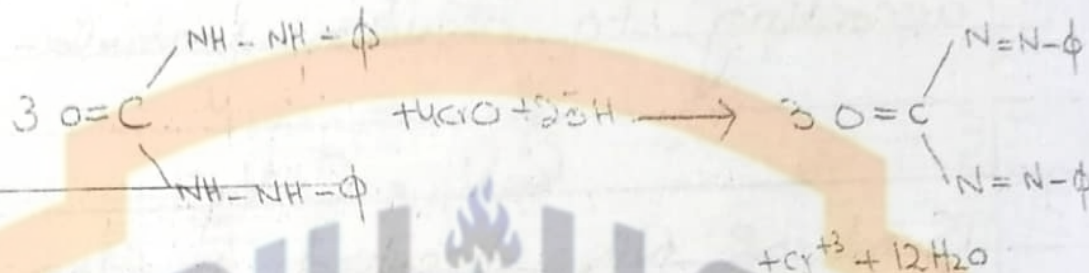

→ $(\text{K}_2\text{Cr}_2\text{O}_7)$

$$\frac{\text{Mol. wt.}}{\text{m. mass}} = \frac{294.18}{52 \times 2} = \frac{294.18}{104} = \frac{2.82}{10} = 0.28 \text{ g/l}$$

→ 0.5% DPC Solution:-

0.25g of DPC in 100ml acetone.

→ 6M H_2SO_4

33ml of H_2SO_4 in 100ml flask.

Experiment No: 28

Determine the amount/L of Cr^{3+} by using Diphenyl carboxide colorimetrically:-

Principle:-

This is colorimetric analysis of Chromium.

Theory:-

Spectroscopy is interaction of radiations with matter. It can be classified on the basis of wavelength and species. Visible spectroscopy is the colorimetry.

Basic requirements of the colorimetric analysis.


Species should be of either soluble in solvent. Suspensions are not to be allowed.

It should be coloured or form the coloured complexes.

Chemistry of Chromium:-


It has atomic no. 24 and group 6th d-block. It belongs to period 4th and transition series. Its atomic wt. is 52 g/mol. Its electronic configuration is,

Material provided by: ece aylin
Arranged by: Umer Abdullah (UA)
GPGC

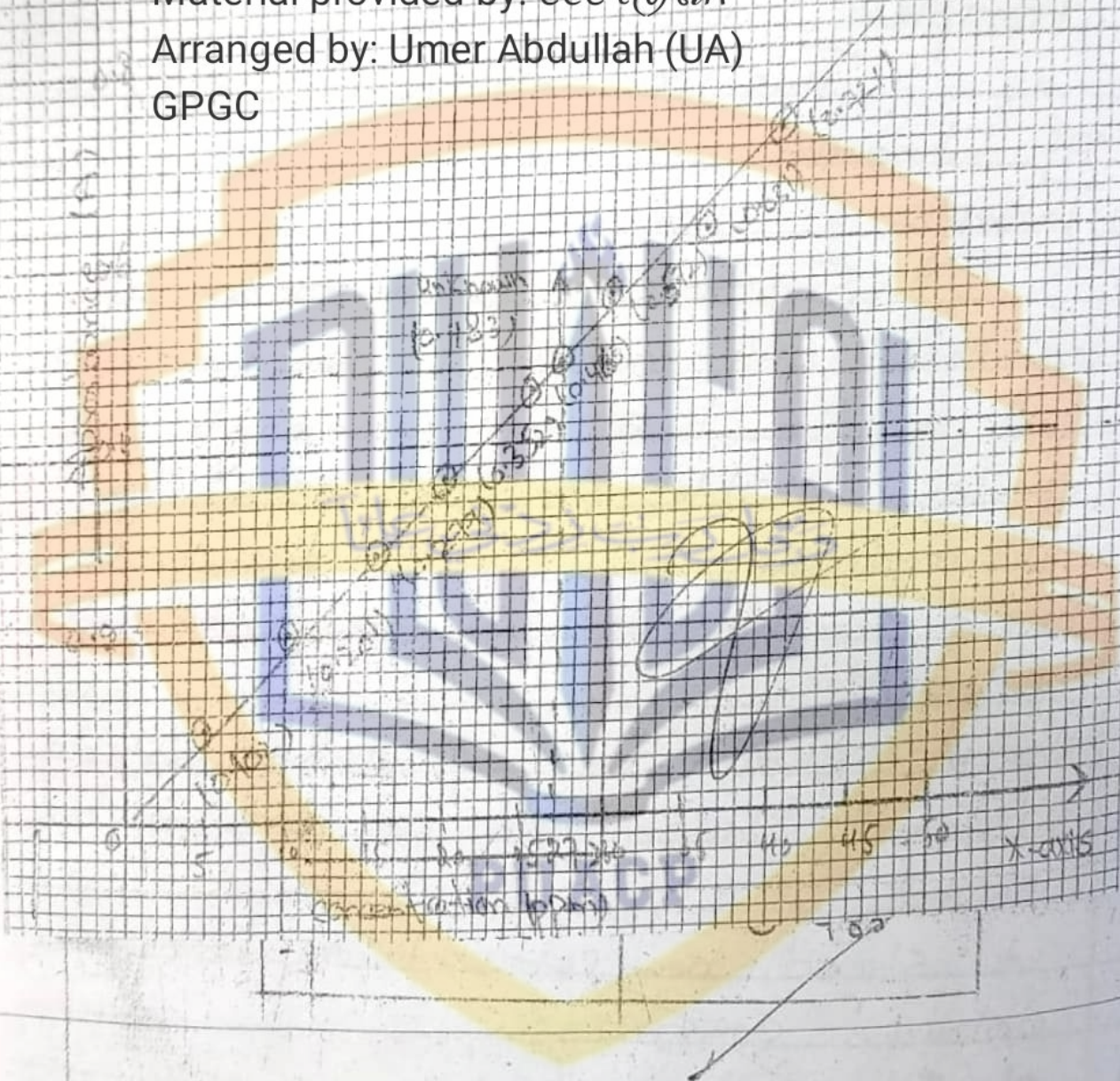

→ Result:-

The conc. of Cr^{+3} in unknown sample is 27.5 ppm.

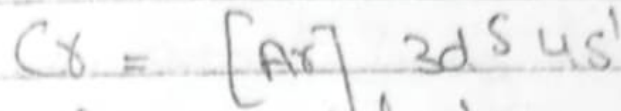

It is steel grey lustrous hard and brittle metal. In larger amount it can be toxic and carcinogenic. The most prominent sample of Cr is the hexavalent Cr^{+6} .

Cr^{III} salts are used for the tanning of leather.

Procedure:-

First of all metal salt soln. of $\text{K}_2\text{Cr}_2\text{O}_7$ in which Cr^{+6} is in $+6$ state was prepared by dissolving 0.28g in 100ml flask and upto mark with distilled H_2O .

0.25% of solution of diphenyl carbazide was prepared by dissolving 0.25g in 100ml acetone and upto mark with distilled H_2O .

Solution was taken in spectronic and reading was taken and graph was plotted b/w absorbance and concentration.