

### Experiment #1

Determine amount per liter of  $\text{Al}^{3+}$   
in given sample gravimetrically.

### Preparation of reagent.

#### Step I:-

Acetic acid

$$M_1 V_1 = M_2 V_2$$

$$0.1 \times 250 = 17 \times V_2$$

$$V_2 = 1.5 \text{ ml}$$

0.1 M Sodium acetate

$$\text{Molarity} = \frac{\text{Molar mass of Solute}}{\text{Molar mass of solvent}} \times \frac{1}{\text{vol.cm}^3}$$

$$0.1 \text{ M} = \frac{x}{82} \times \frac{1}{1000}$$

$$x = 8200$$

**PUACP**

65 ml of acetic acid + 35 ml of sodium acetate =

Buffer solution (pH = 4-5)

#### Step II

2% of 8-hydroxyquinoline

2 g of 8-hydroxyquinoline in 10 ml of ethanol  
and 88 ml of 1 molar acetic acid in dilute soln.

# Experiment No. 1

Determine the amount poss  
of  $\text{Al}^{3+}$  in given Sample gravimetrically.

Theory:-

## Basic Principle


It is a microscopic analytical technique. The word gravimetry is derived from two Greek words Grav means weight and metry means measurements. In this technique, the amount of soluble species is determined by converting it into insoluble forms by using some precipitating agents.

## Use of gravimetric technique

- This technique may be used in wide scale in fertilizing industry for determination of nutrients like  $\text{Zn}$ ,  $\text{Pb}$ ,  $\text{Co}$  etc.
- In cement industry for determination of silica, aluminium and calcium.

In electroplating industry for determination of  $\text{Zn}^{2+}$ ,  $\text{Co}^{2+}$ ,  $\text{Ag}^+$  etc.

## Chemical Equation:-

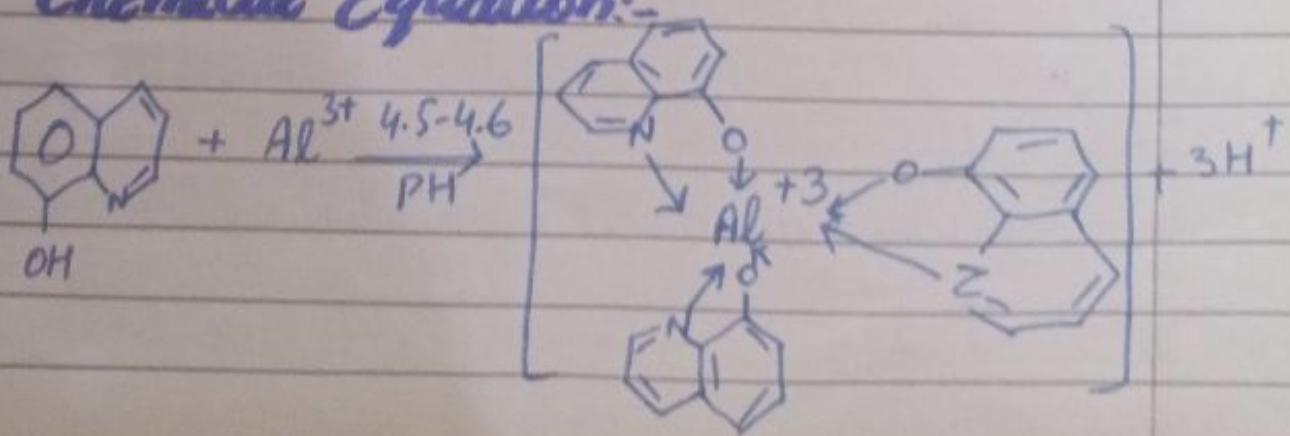

Step III :-

$$1 \text{ molar acetic acid} = V_1 = 100 \text{ cm}^3$$

$$M_1 V_1 = M_2 V_2$$

$$1 \times 100 = 17 V_2$$

$$V_2 = 5.88 \text{ ml}$$

Calculations :-

$$\text{Weight of filter paper} = w_1 = 0.99$$

$$\text{wt. of filter paper + ppt} = w_2 = 1.12$$

$$\text{Net weight} = w = 0.22 \text{ g}$$

$$459 \text{ g of complex contains } Al^{3+} = 27 \text{ g}$$

$$1 \text{ g of complex contains } Al^{3+} = \frac{27}{459}$$

$$0.22 \text{ g of complex contains } Al^{3+} = \frac{27}{459} \times 0.22 = 0.01$$

$$20 \text{ ml of complex contains } Al^{3+} = 0.01 \text{ g}$$

$$1 \text{ ml of complex contains } Al^{3+} = \frac{0.01}{20}$$

$$1000 \text{ ml of complex contains } Al^{3+} = \frac{0.01}{20} \times 1000$$

$$= 0.5 \text{ g/l}$$

Result :-

The given sample contain 0.5 g/l  
of  $Al^{3+}$

## Preparation of reagents:-

For Buffer:- Acetic acid (0.1M)

Step: 1

Given : Required

$$M_1 V_1 = M_2 V_2$$

$$0.1 \times M \times 250 = 17M \times V_2$$

$$1.5 \text{ ml} = V_2$$

Sodium acetate (0.1M)

Molarity =  $\frac{\text{mass of solute}}{\text{molar mass solvent}} \times \text{volume in dm}^3$

$$0.1 = \frac{\gamma}{82} \times \frac{1}{1000}$$

$$\gamma = 8200$$

65ml of acetic acid + 35ml of sodium acetate  
= Buffer soln (PH = 4.5)

Step 2:-

PUACP

2% of  $\beta$  hydroxy quinoline

2g of  $\beta$ -hydroxy quinoline in 10 ml ethanol

and 88 ml of 1 M acetic acid

$$1 \text{ M acetic acid} = V_1 = 100 \text{ cm}^3$$

$$M_1 V_1 = M_2 V_2$$

$$1 \times 1000 = 17 \times V_2$$

$$V_2 = \frac{1000}{17}$$

$$V_2 = 58.8 \text{ ml}$$

## Procedure:-

- Take 20 ml of saturated soln of  $\text{Al}(\text{SO}_4)_2$  in  $\text{H}_2\text{O}$  at  $25^\circ\text{C}$ .
- Add few drops of buffer solution.
- Now add 3-4 drops of buffer solution 8-hydroxy quinoline soln drop wise in above soln.
- Heat it so that precipitates are coagulated.
- Now filter it and weight of filter paper is 0.2 g.

## Use of aluminium

- Aluminium is malleable and ductile metal.
- Alloy of Al is very important. In construction of aeroplanes and other transports.
- Al is used in huge variety of products including foils, lotion, kitchen utensils and window.
- It is good conductor of heat and light.
- It can be polished to give highly refractive surface.

## Experiment # 2

Determine amount/l of  $Zn^{2+}$   
in a given sample gravimetrically.

### Preparation of reagent:-

- 3% solution of anthranilic acid
- 5g of  $NaHCO_3$  were dissolved in 100ml  $H_2O$  to this soln and 3g of anthranilic acid
- 1% metal salt soln, 1g of  $ZnSO_4 \cdot 7H_2O$  dissolved in 99 ml of  $H_2O$

373.39 wt. of complex contain  $Zn^{2+}$  = 65.37 g

$$1 \text{ g} \quad " \quad " \quad " \quad " = \frac{65.37}{373.39}$$

$$0.263 \text{ g} \quad " \quad " \quad " \quad " = \frac{65.37}{373.39} \times 0.26$$

20 ml of sample contain = 0.05

$$1 \text{ ml} \quad " \quad " \quad " \quad " = \frac{0.05}{20}$$

$$1000 \text{ ml} \quad " \quad " \quad " \quad " = \frac{0.05}{20} \times 1000$$

$$= 2.5 \text{ g/l}$$

The given soln contain  $Zn^{2+}$  2.5 g/litre

## Experiment No. 2

Determine amount/liter  
of  $Zn^{2+}$  in a given sample by  
gravimetry.

### Principles-


It is a microscopic analytical technique. The word gravimetry derived from Greek word grav means weight and metry means measurements.

In this technique amount of double species is determined by converting it into insoluble form by using precipitating agents which may be organic or inorganic in nature.

### Use of gravimetric technique

- It is used on ~~wide scale~~ in fertilizer industry.
- For determination of nutrients.
- In cement industry for the determination of silica,  $Al^{3+}$ ,  $CO^{2+}$  etc.
- In electroplating industry for determination of  $Zn^{2+}$ ,  $Co^{2+}$ ,  $Ag^+$  etc.

## Chemical Equation:-

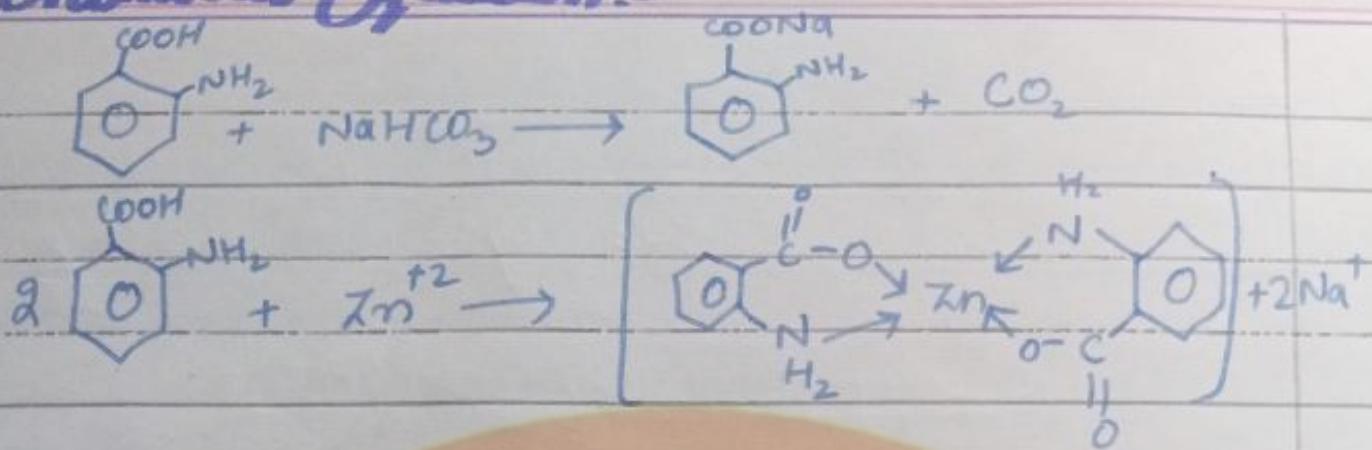

## Preparation of Reagents:-

8% soln of anthracic acid and 5g of  $\text{NaHCO}_3$  were dissolved in 100ml at  $25^\circ\text{C}$  of  $\text{H}_2\text{O}$ . In this soln, add 3g of anthracic acid.

1% Tribasic salt soln of  $\text{ZnSO}_4 \cdot 7\text{H}_2\text{O}$  dissolved in 99ml of water.

### Procedure:-

- Take 20ml of  $\text{ZnSO}_4 \cdot 7\text{H}_2\text{O}$  at  $25^\circ\text{C}$ .
- Now add anthracic acid and 10ml of  $\text{H}_2\text{O}$ .

~~PUACP~~

- Heat it so that ppt are coagulated.
- Now filter it on filter paper and weighed it.
- Weight of filter paper is 0.93 g.
- ppt formed and white in color.

### Uses of Zinc.

- It is used for production of zinc base alloys such as brass, nickel, silica, aluminium oxide.

### Experiment #3

Determine amount per litre of  $\text{Cd}^{2+}$  in given sample by gravimetrically.

### preparation of reagent

- 3% solution of anthranilic acid
- 5g of  $\text{NaHCO}_3$  were dissolved in 100 ml of  $\text{H}_2\text{O}$
- To this soln add 3g of anthranilic acid.
- 1% metal salt i.e. 1 g of  $\text{CdSO}_4$  dissolve in 99 ml of  $\text{H}_2\text{O}$ .

### Calculation:-

$$\text{wt. of filter paper} - w_1 = 1.21$$

$$\text{wt. of filter paper ppt.} - w_2 = 1.43$$

$$\text{Net weight} = 1.43 - 1.21$$

$$= 0.22 \text{ g}$$

$$384 \text{ g of complex contain } \text{Cd}^{2+} = 108 \text{ g}$$

$$1 \text{ g. } \therefore \therefore \therefore = \frac{108}{384}$$

$$0.22 \text{ g of " " " " } = \frac{108}{384} \times 0.22$$

## Experiment No. 3

Determine amount  $\text{ll}$  of  $\text{Cd}^{2+}$  in a given Sample by gravimetry

### Principle:-

In this technique the amount of soluble species is determined by converting it into insoluble species by adding organic or inorganic reagent.


The amount of soluble species is calculated from stoichiometric calculations using knowledge of molar weight of ppt. obtain.

### Use of gravimetric technique

It is used on wide scale on fertilizer industry for determination of nutrients like  $\text{Zn}^{2+}$ ,  $\text{Cd}^{2+}$ ,  $\text{Pd}^{2+}$ ,  $\text{Ag}^+$ .

In cement industry for determination of silica,  $\text{Al}^{3+}$ ,  $\text{Ca}^{2+}$ .

### Chemical Equation

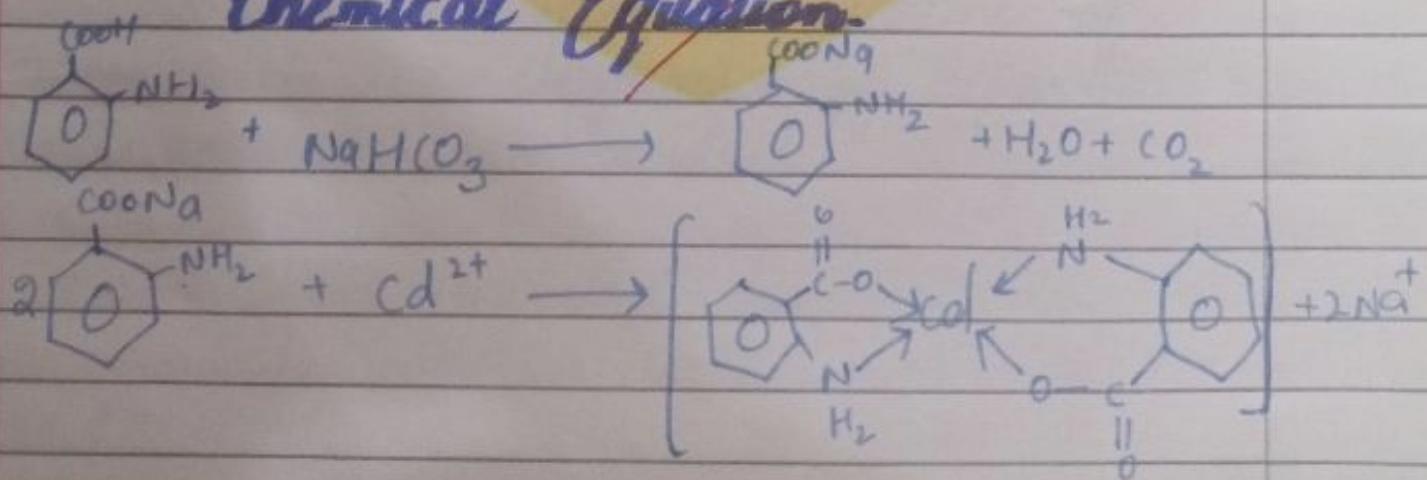

20 ml of sample contain  $\text{cd}^{2+} = 0.669$

$$1 \text{ ml } " " " = \frac{0.06}{20}$$

100ml of sample contain  $\text{cd}^{2+}$

$$= \frac{0.06}{20} \times 1000$$
$$= 3.02 \text{ g/l}$$

Result:-

The given sample contain  
 $\text{cd}^{2+} 3.02 \text{ g/l}$  علی علی ذریعہ ذرا سب

PUACP

# Preparation of Organic Reagent

8% of anthracic acid in  
 $\text{NaHCO}_3$ , 5g of  $\text{NaHCO}_3$  dissolve in 100ml  
of  $\text{H}_2\text{O}$  and add 3g of anthranilic acid.

1% of metal soln, 1g of  $\text{Cd}(\text{NO}_3)_2$   
dissolve in 100ml of  $\text{H}_2\text{O}$ .

## Procedure:-

Take 20 ml calculated solution of  
 $\text{Cd}(\text{NO}_3)_2$  at 25°C.

Now add anthranilic acid and add  
10 ml of  $\text{H}_2\text{O}$ .

Heat it so that precipitates are  
coagulated.

Now filter it on filter paper,  
light pink ppt were obtained.

## Uses of $\text{Cd}^{2+}$

The metal is ~~plastic~~ component of  
alloys with low-coefficient of expansion  
and good fatigue resistant.

$\text{Cd}^{2+}$  is important for batteries.

It is used in some control rods  
and shield in nuclear reaction.

## Experiment #4

Determine amount / g of  
 $\text{Co}^{2+}$  in given sample gravimetrically.

### Preparation of Reagent :-

- 3% soln of anthranilic acid
- 5g of  $\text{NaHCO}_3$  of anthranilic acid
- 1 g of metal salt + 1g of cobalt nitrate in 99 ml  $\text{H}_2\text{O}$ .

### Calculations :-

$$\begin{aligned}\text{wt. of filter paper} &= w_1 = 0.93 \\ \text{wt. of filter paper} + \text{ppt} &= w_2 = 1.61 \\ \text{net weight} &= w = 1.61 - 0.93 \\ &= 0.23\end{aligned}$$

331 g of metal complexes contain  $\text{Co}^{2+} = 59$

$$1 \text{ g } " = = = \frac{59}{331}$$

$$0.23 \text{ g } " = = = \frac{59}{331} \times 0.23 \\ = 0.049$$

# Experiment # 4

Determine amount  $\text{M}$  of  $\text{Co}^{2+}$  in given gravimetry

## Principle:-

In this technique the amount of soluble compound is determined by using some organic or inorganic reagents.


The amount of soluble specie is calculated from stoichiometric calculations using knowledge of precipitates obtained.

## Uses of gravimetric technique

It is used on wide scale in fertilizer industry for determination of nutrients like  $\text{Zn}^{2+}$ ,  $\text{Cd}^{2+}$ ,  $\text{SO}_4^{2-}$ ,  $\text{Pd}^{2+}$ .

In Cement industry for determination of Silica, alumina etc.

## Chemical Equation

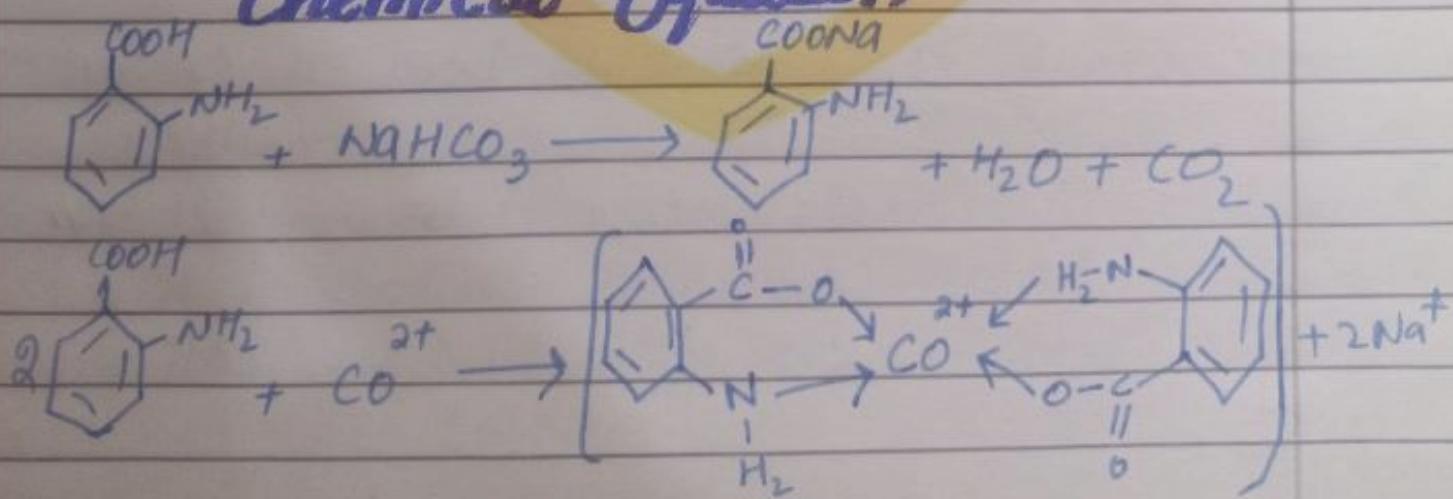

20 ml of sample contain  $\text{Co}^{2+} = 0.04$

$$1 \text{ ml} \quad " \quad " \quad " = \frac{0.04}{20}$$

$$1000 \text{ ml} \quad " \quad " \quad " = \frac{0.04}{20} \times 1000 \\ = 2.05$$

Result :-

The given sample contain  $\text{Co}^{2+}$   
is 2.05 g/l.

جامعة زنجبار

PUACP

## Preparation of Reagents

- 5g of  $\text{NaHCO}_3$  dissolve in 100ml  $\text{H}_2\text{O}$  and add 8g of anthranilic acid.
- 1% metal salt solution 1g of cobalt nitrate dissolve in 999 ml of  $\text{H}_2\text{O}$ .

## Procedure:-

- Take 20ml of soln of calculated cobalt nitrate at 25°C.
- Now add anthranilic acid till ppt obtained.
- Heat ~~soln~~ to coagulate ppt.
- cool solution at room temperature, weight filter paper and from it calculate weight of ppt formed.

## Uses of $\text{Co}^{3+}$ P.D.A.C.P

- It is used as a source of gamma radiation.
- It is used in razor blades.
- It is used as permanent magnet.
- It is used as alloy for engines and gas turbines.

## Experiment # 5:-

Determine amount/l of  
Ni<sup>2+</sup> gravimetrically.

### Preparation of reagents:-

⇒ 1% NiSO<sub>4</sub> soln

and 1 g of NiSO<sub>4</sub> is dissolved in H<sub>2</sub>O  
make volume upto 100 ml.

⇒ 1% Salicylaldoxime

Dissolve 1 g of Salicylaldoxime in 10-15 ml  
of ethanol for complete solubility and  
then dilute the soln upto 100 ml in H<sub>2</sub>O.

### Calculations

**PUACP**

$$\text{wt. of filter paper} = w_1 = 1.19$$

$$\text{wt. of filter paper + ppt} = w_2 = 1.34$$

$$\begin{aligned}\text{Net weight} &= w = 1.34 - 1.19 \\ &= 0.15 \text{ g}\end{aligned}$$

## Experiment No. 5

Determine amount/l.  
of  $\text{Ni}^{2+}$  gravimetrically.

Principle:-

Theory:- It is microscopic analytical technique.


In this technique amount of soluble species is determined by converting it into insoluble form using organic or inorganic precipitating agent.

Uses of gravimetric technique

It is used on wide scale for determination of nutrients like  $\text{Zn}^{2+}$ ,  $\text{Co}^{2+}$ ,  $\text{Pb}^{2+}$ ,  $\text{SO}_4^{2-}$ .

In cement industry for determination of silica,  $\text{Al}^{3+}$ ,  $\text{Co}^{2+}$ . In electroplating industry for determination of  $\text{Zn}^{2+}$ ,  $\text{Co}^{2+}$ ,  $\text{Ag}^{2+}$ .

Chemical Equation

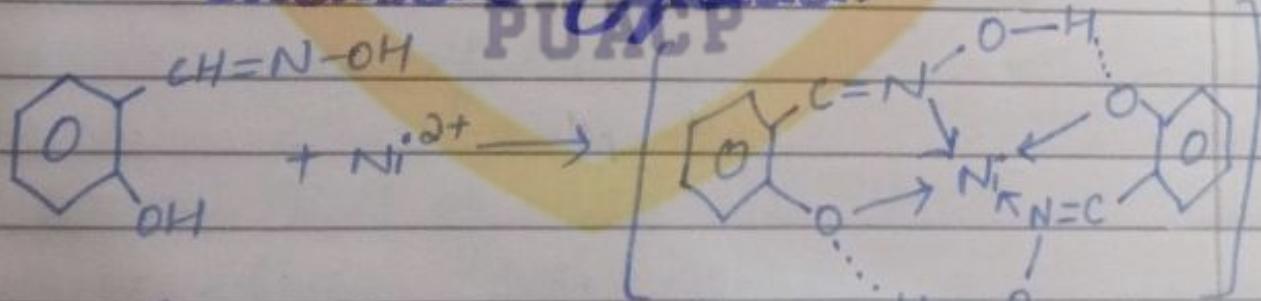

Preparation of Reagents

1%  $\text{NiSO}_4$  soln:-

$$330 \text{ g of complex contain } Ni^{2+} = 58.6 \text{ g}$$

$$\frac{1 \text{ g}}{330 \text{ g}} = \frac{58.6}{x} \Rightarrow x = \frac{58.6}{330} \times 0.15$$

$$= 0.02 \text{ g}$$

$$20 \text{ ml of } \sim Ni^{2+} = \frac{0.02}{20}$$

$$\frac{1 \text{ ml}}{1000 \text{ ml}} = \frac{0.02}{20} \times 1000$$

$$= 1.33 \text{ g/l}$$

Result :-

Sample The amount of  $Ni^{2+}$  in given  
is 1.33 g/l. دلتا

PUACP

1 g of  $\text{NiSO}_4 \cdot 6\text{H}_2\text{O}$  is dissolved in 99 ml of water and make volume upto 100 ml.

### Salicyl aldoxime :-

- Dissolve 1 g of salicylaldoxime ion in 10-15 ml of ethanol for complete solubility and dilute soln upto 100 ml of  $\text{H}_2\text{O}$ .

### Procedure:-

- Add 20 ml of  $\text{NiSO}_4$  in a flask. Add few drops of salicylaldoxime soln until ppt are formed.
- Heat soln that ppt are coagulated.
- Green color ppt are formed.
- filter out and weight ppt obtained.

### uses of $\text{Ni}^{2+}$

- Most uses of nickel are alloys metals namely nickel salt, bronze and nickel iron.
- stainless steel is made of nickel.

## Experiment #6 :-

Determine amount of  $\text{Ni}^{2+}$  in given sample using nitron gravimetrically.

### Preparation of organic Reagent:-

$\Rightarrow$  0.1 % soln of nitron

For 0.1 % soln of nitron takes 0.1 g of nitron and dissolve in 5M acetic acid and make total soln 100ml.

$\Rightarrow$  5 ml Acetic-acid

Take 3.5 ml of acetic acid. Dissolve in 95 ml of water to make 5M acetic acid.

### Calculation :-

**PUACP**

$$\text{wt. of filter paper} = w_1 = 1.5 \text{ g}$$

$$\text{wt. of filter paper + ppt} = w_2 = 1.71$$

$$\begin{aligned}\text{Net weight} &= w = w_2 - w_1 = 1.71 - 1.5 \\ &= 0.21\end{aligned}$$

## Experiment No. 6

Determine amount of  $\text{NO}_3^-$  in given sample using nitron by gravimetrically.

### Principle:-

Gravimetry is microscopic analytical technique.

### Theory:-

In this technique the amount of soluble species is determined by converting into insoluble species using organic or inorganic reagent.

### Preparation of solution

For 10% soln of nitron, take 0.1g of nitron and dissolve in 5M of acetic acid and make total volume of soln 100ml.

If nitron soln does not dissolve readily we can also dissolve few ml of ethanol in 95 ml total volume of 100ml.

### 5M acetic acid

Take 3.5M of acetic acid and dissolve in 95 ml of  $\text{H}_2\text{O}$  to make 5M acetic acid.

### 0.1% solution of $\text{Ag}(\text{NO}_3)_2$

In distilled water dissolve 0.1g of  $\text{Ag}(\text{NO}_3)_2$

374 g of complex contain  $\text{NO}_3^- = 62 \text{ g}$

$$1 \text{ g } " - " = \frac{62}{374}$$

$$0.21 \text{ g } " - " = \frac{62}{374} \times 0.21$$

20 ml Sample contain  $\text{NO}_3^- = 0.034$

$$1 \text{ ml } " - " = \frac{0.034}{20}$$

$$1000 \text{ ml } " - " = \frac{0.034 \times 1000}{20} = 1.75 \text{ g/l}$$

Result:-

The amount of  $\text{NO}_3^-$  in given sample is 1.75 g/l.

in 5g of  $\text{CH}_3\text{COOH}$  to make total 100ml of soln.

## Chemicals

- Ethanol
- acetic acid
- Nitron
- Aluminium nitrate


## Uses of nitrates:-

A rich source of nitrates in human body comes from leafy food.

It is used in fertilizer and agriculture industry.

A rich source is found in Earth crust and in the form of ore.

## Equation:-

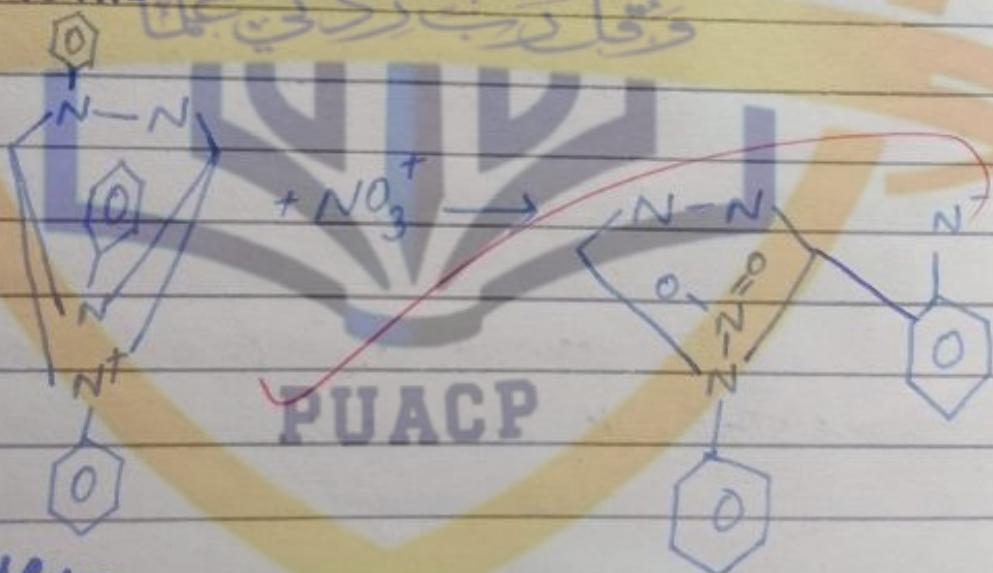

## Procedure:-

- Take 20 ml of saturated of given sample at 25°C.
- Add nitron soln and few drops of ethanol to dissolve the nitron completely.
- Heat it so that ppt are being formed and coagulate filter and weight it.

## Experiment # 7 :-

Determine amount /g of  $\text{Bi}^{3+}$  gravimetrically.

⇒ Preparation of reagent:-

Dissolve 1 g of  $\text{Bi}(\text{NO}_3)_3$  in  $\text{H}_2\text{O}$  and make it upto 100 ml.

⇒ 1% Pyrogallol

Dissolve 1 g of Pyrogallol in hot  $\text{H}_2\text{O}$  and dilute upto 100 ml.  $\text{H}_2\text{O}$  used because Pyrogallol is insoluble in cold  $\text{H}_2\text{O}$ .

Calculations :-

$$\text{wt. of filter paper} = w_1 = 1.27$$

$$\text{wt. of filter paper + ppt} = w_2 = 1.37$$

$$\text{Net weight} = w = w_2 - w_1 = 1.37 - 1.27$$

$$331 \text{ g of complex contain } \text{Bi}^{3+} = 0.1 \text{ g}$$

$$1 \text{ g " " " } = \frac{209}{331}$$

$$0.1 \text{ g " " " } = \frac{209}{331} \times 0.1$$

# Experiment No. 7

Determine amount of  $\text{Bi}^{+3}$  gravimetrically.

## Principle:-

It is microscopic analytical technique

## Theory:-


In this technique the amount of soluble species is determined by converting it into insoluble species using organic or inorganic reagents.

## Chemical:-

1%  $\text{Bi}(\text{NO}_3)_3$  soln

1% Pyrogallol soln

## Equation:-

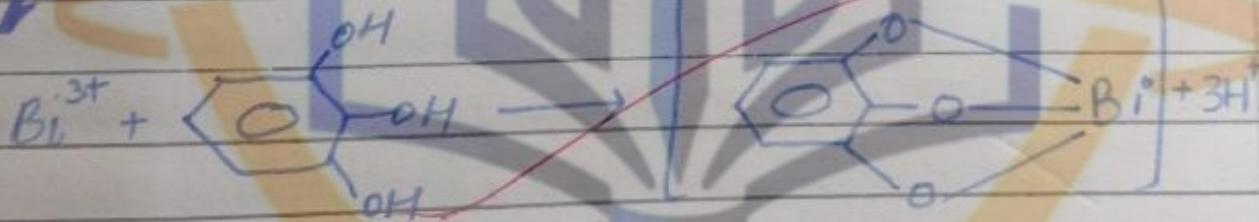

## Preparation of soln:-

1%  $\text{Bi}(\text{NO}_3)_3$

Dissolve 1g of  $\text{Bi}(\text{NO}_3)_3$  in water and make it dilute upto 100 ml.

1% Pyrogallol

Dissolve 1g of pyrogallol in hot  $\text{H}_2\text{O}$  and make volume upto 100ml. Pyrogallol is insoluble in cold water.

20 ml sample contain  $\text{Bi}^{3+} = 0.06$


$$1 \text{ ml} \quad " \quad " \quad " = \frac{0.06}{20}$$

$$1000 \text{ ml} \quad " \quad " \quad " = \frac{0.06}{20} \times 1000$$

$$= 3.15 \text{ g/l}$$

Result :-

The amount of  $\text{Bi}^{3+}$  in a given  
sample is 3.15 g/l.

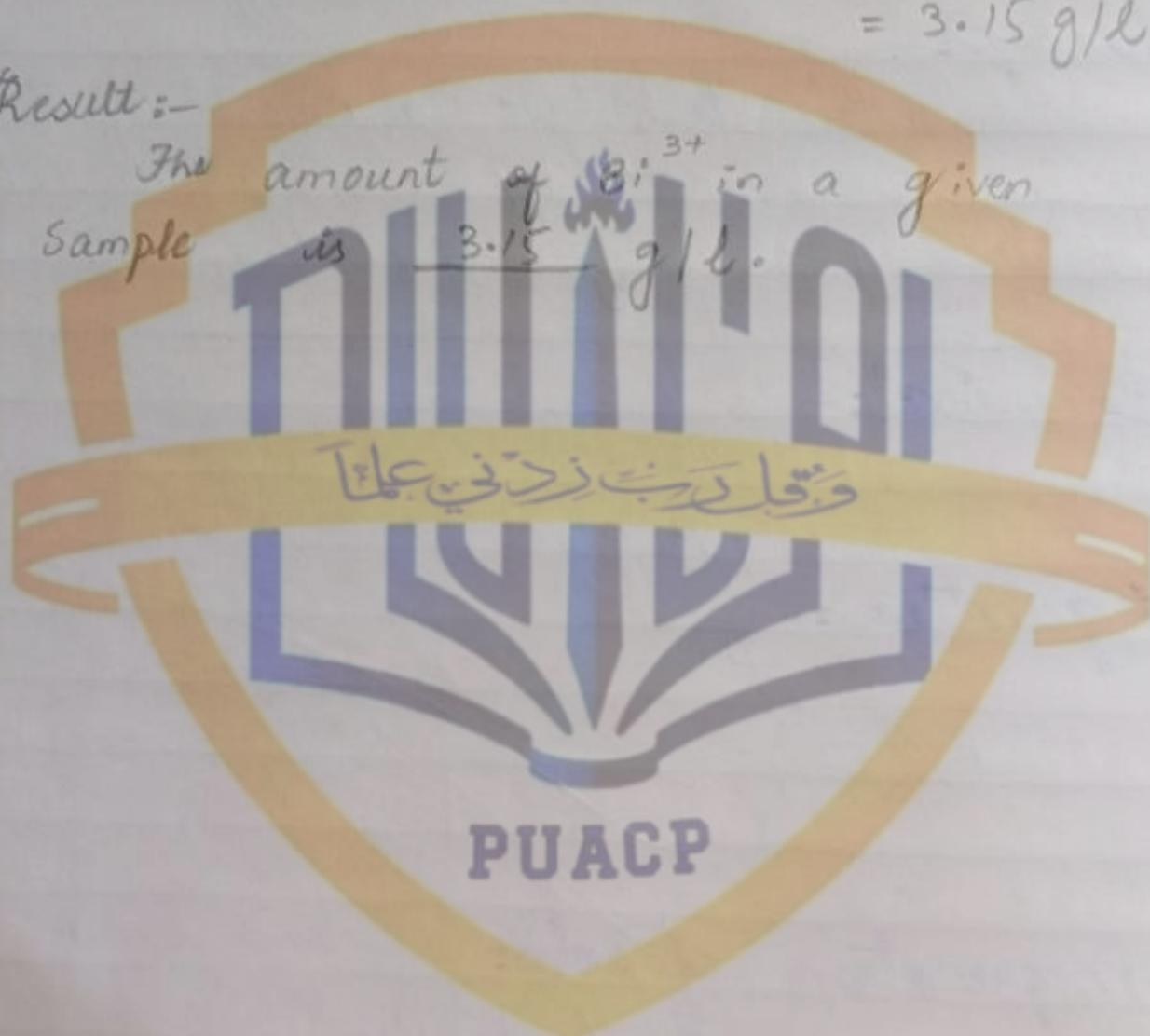

## Uses of Bi<sup>3+</sup>

- It is used as thermocouple materials.
- It is used in cosmetics.
- The alloy of Bismuth is used in auto-matic fine sprinkles.
- In detection system and in electronic fuses.
- Bi<sup>3+</sup> is usually mixed with other metal like lead, iron or cd to form melting alloy.

## Procedure:-

Take 20 ml of sample solution. Add 1% Pyrogallol soln dropwise in the sample soln until in soln ppt are formed.

Heat the soln so that ppt are coagulated.

No filter it and weight ppt formed.

**PUACP**

Brown color ppt are formed.

## Experiment # 8 :-

Determine amount of  $\text{Fe}^{3+}$  in given sample gravimetrically.

Preparation of Reagent :-

- Buffer soln:-

$$M_1V_1 = M_2V_2$$

$$0.1 \times 250 = 17 \times V_2$$

$$V_2 = 1.5 \text{ ml}$$

- 0.1 M sodium Acetate :-

65 ml of  $\text{CH}_3\text{COOH} + 35 \text{ ml of } \text{CH}_3\text{COONa} =$

Buffer Solution (PH 4-5)

- 2% 8 hydroxyquinoline

2 g of 8 hydroxyquinoline dissolve in 10ml ethanol for complete solubility and 88 ml of 1 M acetic acid is added.

1 M  $\text{CH}_3\text{COOH}$

$$M_1V_1 = M_2V_2$$

$$1 \times 100 = 17 \times V_2$$

$$V_2 = 5.8$$

# Experiment # 8

Determine amount/l of  $\text{Fe}^{3+}$  in a given sample gravimetrically.


Principle:-

It is microscopic analytical technique.

Theory:-

In this technique the amount of sample species is determined by converting it into insoluble species using some organic or inorganic reagents.

Chemical equation:-

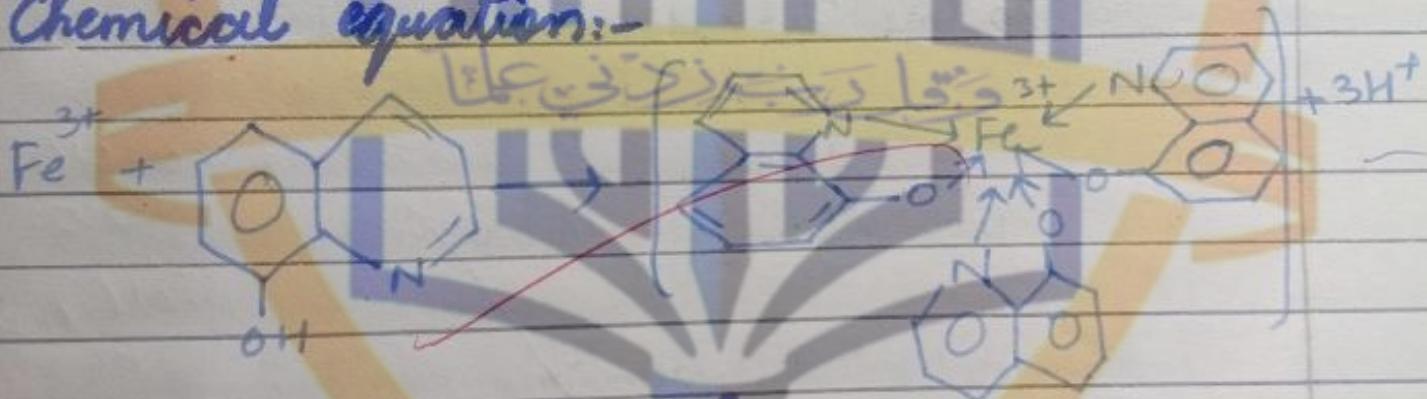

Preparation of soln

Buffer Soln  
Given = Required

$$M_1 V_1 = M_2 V_2$$

$$0.1 \times 250 = 17 \times V_2$$

$$V_2 = \frac{0.1 \times 250}{17}$$

$$V_2 = 1.5 \text{ ml}$$

### Calculation:-

$$\begin{aligned} \text{wt. of filter paper} &= w_1 = 1.20 \text{ g} \\ \text{wt. of filter paper + ppt} &= w_2 = 1.34 \text{ g} \\ \text{Net weight} &= w = 1.34 - 1.20 \\ &= 0.14 \text{ g} \end{aligned}$$

488 g of complex contain  $\text{Fe}^{3+} = 55.8$

$$\begin{aligned} 0.14 \text{ g} &= \frac{55.8}{488} \\ &= \frac{55.8}{488} \times 0.14 \end{aligned}$$

20 ml of complex contain  $\text{Fe}^{3+} = 0.016$

$$\begin{aligned} 1000 \text{ ml} &= \frac{0.016}{20} \\ &= \frac{0.016}{20} \times 1000 \end{aligned}$$

**PUACP** = 0.8 g/l

The amount of  $\text{Fe}^{3+}$  in a given sample is 0.8 g/l.

### Result:-

0.1 M  $\text{CH}_3\text{COONa}$

Molarity =  $\frac{\text{mass of solute}}{\text{molar mass}} \times \frac{1}{\text{vol dm}^3}$

$$0.1 = \frac{x}{82} \times \frac{1100}{1000}$$

$$x = 0.1 \times 82 \times 0.1$$

$$x = 0.82 \text{ g}$$

65 ml of  $\text{CH}_3\text{COOH}$  + 35 ml of  $\text{CH}_3\text{COONa}$  = Buffer soln  
PH=4.5

2% 8 hydroxy quinoline

~~PUACP~~

2g of hydroxy quinoline dissolved in 100 ml  
of ethanol, for example solubility 88 ml of  
1M acetic acid is added

1 M acetic acid

Given = Required

$$M_1 V_1 = M_2 V_2$$

$$V_2 = 5.88 \text{ ml}$$

## Procedure:-

- Take 20 ml of saturated solution of  $\text{Fe}^{3+}$  at 25 °C.
- Add few drops of buffer solution.
- Now add 8-hydroxyquinoline soln dropwise in above soln.
- Heat it so that ppt are coagulated.
- Green ppt are obtained.

دُقَارَبَ زَرْنَقِي عَلَى

PUACP

### Experiment # 9:-

Determine amount/l. of  $\text{Cl}^-$  using  
0.05 M  $\text{AgNO}_3$  soln.

$\Rightarrow$  Preparation of reagent:-

(a) 0.02 M of  $\text{NaCl}$  :-

$$0.025 \text{ M} = \frac{x}{59} \times \frac{1}{100/1000}$$

$$x = 0.025 \times 59$$

$$x = 1.478 \text{ g/100 ml}$$


(b) 0.05 M of  $\text{AgNO}_3$

$$0.05 \text{ M} = \frac{x}{169} \times \frac{1}{100/1000}$$

$$x = \underline{\underline{0.84}}$$

$$x = 0.84 \text{ g/100 ml}$$

Chemical Equation :-

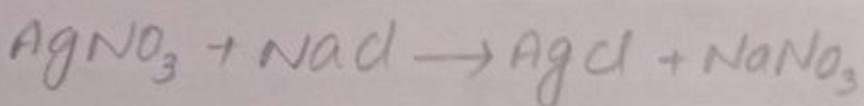

# Experiment No. 9

Determine amount/l of  $\text{Cl}^-$  using 0.05 M  $\text{AgNO}_3$  soln.

## Principle:-

It is precipitation titration.

## Theory:-

The reciprocal of resistance is termed as conductance.

Unit of conductance is Siemen. (S).

Resistance is directly proportional to length and inversely proportional to area.

$$R = \rho \frac{l}{A}$$

$\rho$  is known as resistivity. Reciprocal of resistivity is called conductivity.


$$K = \frac{1}{\rho} \text{PUACP}$$

Conductance of electrolyte soln depends upon ions present and their concentration at particular temperature.

When soln is diluted, conductance decreases since few ions are present. For strong electrolytes molar conductivity increases.

vol. of AgNO <sub>3</sub>	conductance (mS)
0	1.47
3	1.314
6	1.190
9	1.083
12	1.172
15	1.266
18	1.36
21	1.516
24	1.621
27	1.718
30	1.80

Calculation:-

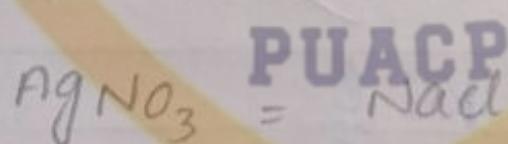

$$M_1 V_1 = M_2 V_2$$

$$0.05 \times 9 = M_2 \times 25$$

$$M_2 = 0.018 \text{ M}$$

Exp No 9.

• Scale

- Along x-axis


$$1 \text{ cm} = 6$$

$$1 \text{ mm} = 0.1$$

- Along y-axis

$$1 \text{ cm} = 1.1$$

$$1 \text{ mm} = 0.11$$

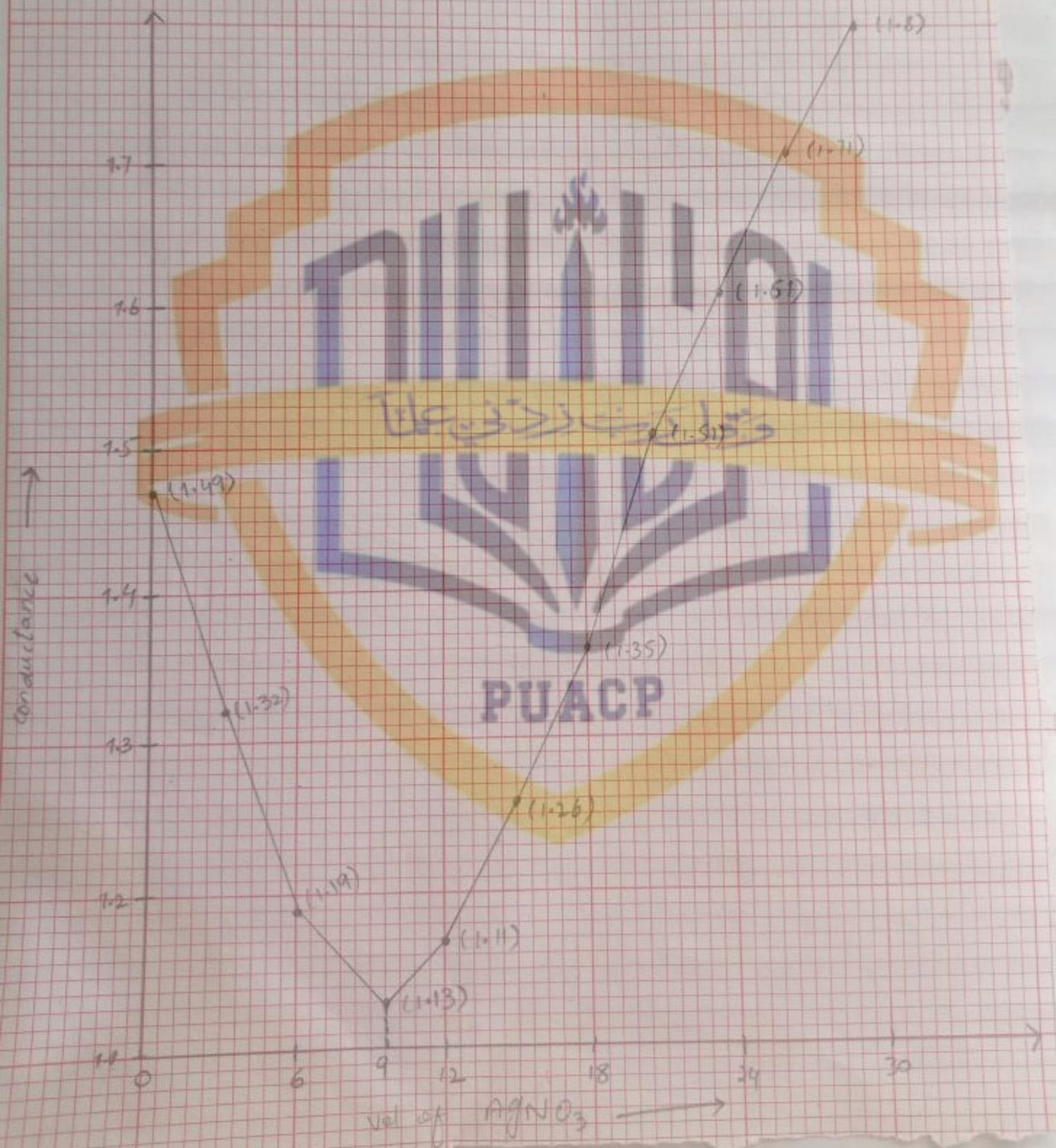

but it approaches to limiting value known as molar conductivity at infinite dilution.

### Conductometer

Instrument that is used to measure conductance. In short conductance decrease due to rapid fall of conc. of  $\text{Cl}^-$  ions. After equivalence point, conductivity rises due to increase in conc. of  $\text{Ag}^+$  ions.

### Uses of $\text{Cl}^-$ ion.

- Chlorine belongs to group VII A of periodic table and have uses.
  - Used for preserving food.
  - used as bleaching agent in paper and pulp industry.

### Uses of $\text{Ag}^+$ ion

- Atomic no. is 47 and mass no is 108.
  - used in jewellery electroplating.
  - In decoration of sweets.
  - Co-ordination no. is 2 and complexes are linear.

### Procedure:-

- 25 ml of sample soln  $\text{NaCl}$  was taken in beaker  $0.05 \text{ M } \text{AgNO}_3$  taken in burette.
- Add 2 ml of it to sample and measure the conductance.

Strength of  $\text{Cl}^-$  = molarity  $\times$  atomic wt. of  $\text{Cl}^-$

$$= 0.018 \times 35.5$$


$$= 0.639 \text{ g/dm}^3$$

Result:-

The amount/l of  $\text{Cl}^-$  in a given sample is  $0.639 \text{ g/l}$ .

PUACP

- The conductance falls in start and rises after equivalence point.
- Plot graph b/w volume of  $\text{AgNO}_3$  on x-axis and conductance on y-axis.

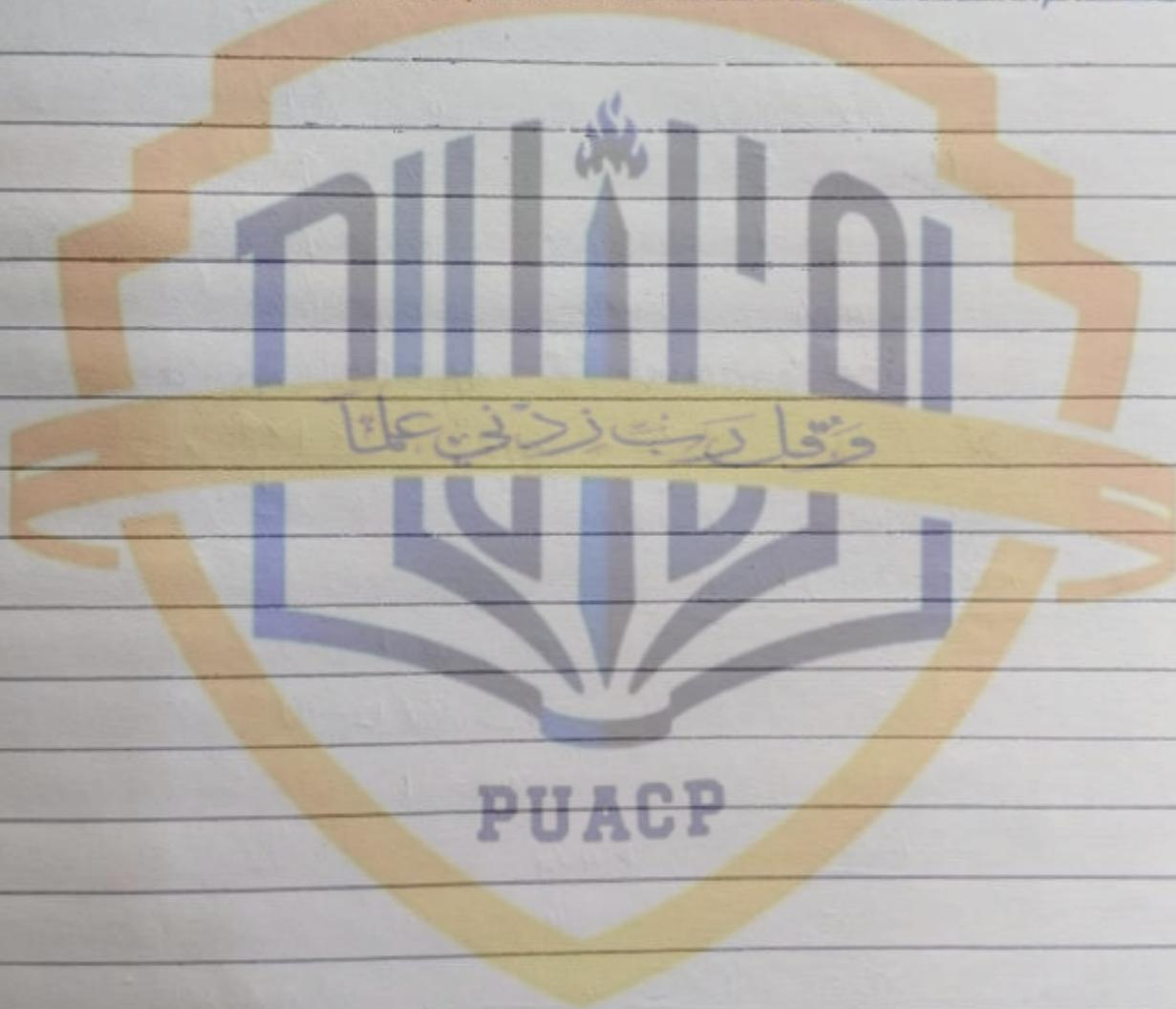

## Experiment # 10:-

Standardize given mixture  
of  $\text{CH}_3\text{COOH} + \text{HCl}$  conductometrically.

### Preparation of reagent:-

(i)



$$M_1 = ?$$

$$V_1 = 5$$

$$n_1 = 1$$

$$\frac{M_1 V_1}{n_1} = \frac{M_2 V_2}{n_2}$$

$$M_2 = 0.1$$


$$V_2 = 9$$

$$n_2 = 1$$

$$M_1 = \frac{M_2 V_2}{n_2} \times \frac{n_1}{V_1} = \frac{0.1 \times 9}{1} \times \frac{1}{5}$$

$$M_1 = 0.18 \text{ M}$$

(ii)

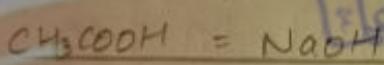

$$M_3 = ?$$

$$V_3 = 5$$

$$n_3 = 1$$

$$M_4 = 0.1$$

$$V_4 = 3(12 - 9)$$

$$n_4 = 1$$

$$\frac{M_3 V_3}{n_3} = \frac{M_4 V_4}{n_4}$$

$$M_3 = \frac{0.1 \times 3}{1} \times \frac{1}{5} = 0.06 \text{ M}$$

(i) Strength of  $\text{HCl}$  = Molarity  $\times$  mol. wt

$$= 0.18 \times 36.5$$

$$= 6.57 \text{ g/dm}^3$$

(ii) Strength of  $\text{CH}_3\text{COOH}$  = Molarity  $\times$  mol. wt

$$= 0.06 \times 60$$

$$= 3.6 \text{ g/dm}^3$$

## Experiment No. 10

Standardize the given mixture  
of  $\text{CH}_3\text{COOH} + \text{HCl}$  conductometrically

Principle:-

It is acid base titration.

Theory:-

The reciprocal of resistance is conductance  
Units of conductance is Siemens (S).

Resistance is directly proportional to  
length and inversely proportional to Area.

$$R = \rho \frac{l}{A}$$


Conductivity is the reciprocal of resistivity  
conductance of electrolyte soln. depend  
upon ion present and their conc. at  
particulate temp.

When soln is diluted  
conductance decreases ~~because~~ few ions are  
present. For strong electrolyte molar conductivity  
increases as dilution increases.

But it approaches limiting value known  
as molar conductivity at infinite dilution.

Basis for conductometric titration

In start, the conductivity decreases  
due to rapid fall of conc. of  $\text{H}^+$  ions of

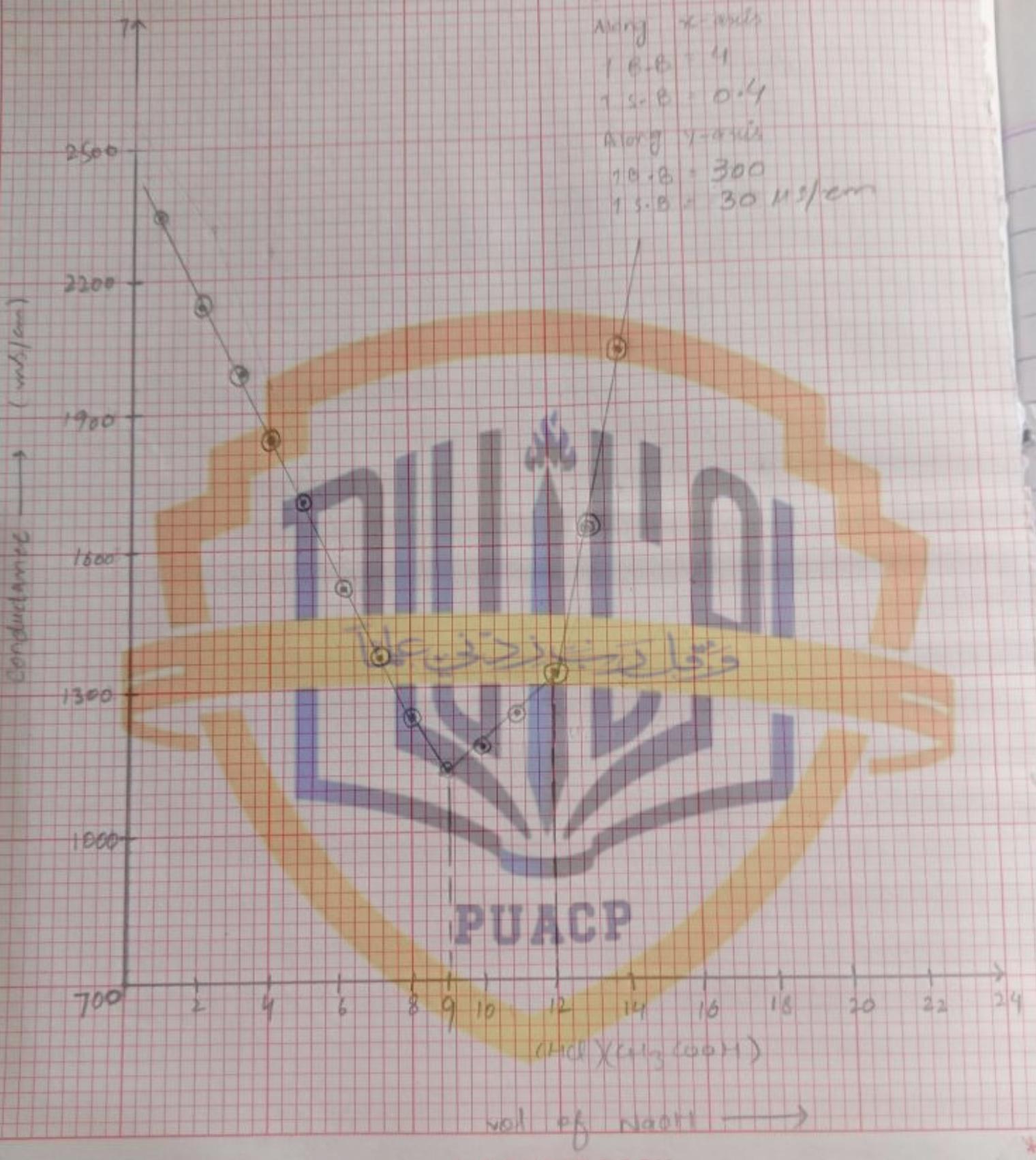

RIZVI NOTE BOOK

(i) HCl	=	6.57 g/dm <sup>3</sup>
(ii) CH <sub>3</sub> COOH	=	3.6 g/dm <sup>3</sup>

## Observation & Calculations:-

	Volume of NaOH	Observed Conductance (μS/cm)
1	1	1990
2	2	1840
3	3	1709
4	4	1545
5	5	1365
6	6	1243
7	7	1190
8	8	1260
9	9	1270
10	10	1320
11	11	1657
12	12	2050
13	13	2230

## Result:-

Strength of each acid in mixture

as

$$(i) \text{ HCl} = 6.57 \text{ g/dm}^3$$

$$(ii) \text{ CH}_3\text{COOH} = 3.6 \text{ g/dm}^3$$

of strong acid HCl, after equivalence the conductivity slowly rises due to slow dissociation of  $H^+$  ions of weak acid. After that  $H^+$  ions of weak acid competes with  $OH^-$  ions of strong base (NaOH) and sharp rise of conductivity is obtained.

### Procedure:-

- Take 25 ml of  $CH_3COOH$  and 25 ml of HCl and takes mixture of both acids in beaker.
- Take 0.1 M NaOH in burette.
- Now 2 ml of it is added to sample soln.
- Measure the conductance. The conductance falls in start and slowly rise after second equivalent point.
- The graph b/w vol. of NaOH and conductance is plotted.

## Experiment # 11:-

Determine conc. of  $\text{Ni}^{2+}$

Calorimetrically or spectrophotometrically:

Preparation of reagent:-

1- Stock soln:- Molar mass of compound:-

Molar mass of  $\text{Ni}^{2+}$

$$= \frac{237.7}{58.7} = 0.4 / 1000 \text{ ml H}_2\text{O}$$

2- Prep of different ppm conc. from stock soln

(i) Prep of 10 ppm soln

Given = Required

$$1000 \times V_1 = 10 \times 1000$$

$$V_1 = 1 \text{ ml / 100 ml H}_2\text{O}$$

(ii) Prep of 15 ppm soln

Given = Required

$$\frac{C_1 V_1}{1000} = \frac{C_2 V_2}{100}$$

$$1000 \times V_1 = 15 \times 100$$

$$V_1 = 1.5 \text{ ml / 100 ml H}_2\text{O}$$

(iii) Prep of 20 ppm soln

$$C_1 V_1 = C_2 V_2$$

$$1000 V_1 = 20 \times 100$$

$$V_1 = 2 \text{ ml / 100 ml H}_2\text{O}$$

## Experiment No. 11

Determine conc. of  $\text{Ni}^{2+}$  calorimetrically or spectrophotometrically.

### Theory:-

Measurement of sample by calibration curve as a result of standard soln analysis is known as calorimetry analysis.

### spectroscopy

It is interaction of electromagnetic radiations with matter.

#### Types:-

##### Atomic : ~~Emission~~

It is further subdivided into atomic absorption spectroscopy.

##### Molecular :-

Further include UV, visible spectroscopy, IR spectroscopy, NMR spectroscopy, mass spectroscopy. PUACP

### Sources of excitation

Sources of excitation can be flamed by hollow cathode lamp. Inductive coupled plasma (ICP).

### Beer Lambert Law

UV/visible spectra based upon beer lambert law which states that

(iv) prep of 25 ppm soln:-

$$C_5V_5 = C_6V_6$$
$$1000 \times V_5 = 25 \times 100$$

$$V_5 = 2.5 / 100 \text{ ml H}_2\text{O}$$

(v) Prep of 30 ppm soln:-

$$C_7V_7 = C_8V_8$$
$$1000 \times V_7 = 30 \times 100$$

$$V_7 = 3 \text{ ml} / 100 \text{ ml}$$

(vi) prep of 35 ppm soln:-

$$C_9V_9 > 5 \rightarrow C_{10}V_{10}$$
$$1000 \times V_9 = 35 \times 100$$

$$V_9 = 3.5 \text{ ml} / 100 \text{ ml}$$

(vii) Prep of 40 ppm soln:-

$$C_{11}V_{11} = C_{12}V_{12}$$

$$1000 \times V_{11} = 40 \times 100$$

$$V_{11} = 4.0 \text{ ml H}_2\text{O}$$

Preparation of DMG

Dissolve 0.1 g of DMG in ethanol and make vol. upto 100ml.


The absorption of light by sample is directly proportional to conc. and path length of sample.

$$A \propto cl$$

$$A = a.c.l$$

As we fix length of path length that is  $l = 1\text{ cm}$  then

$$A \propto c$$

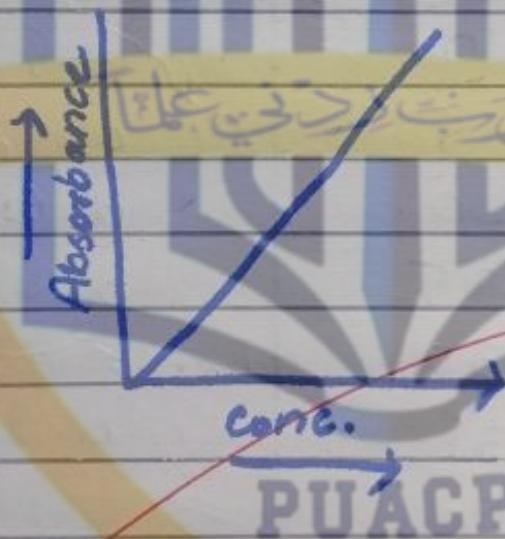


Source of light:-

Deutrium tungsten.

Cuvette

Cuvettes are made up of quartz

Chemical Equation:-

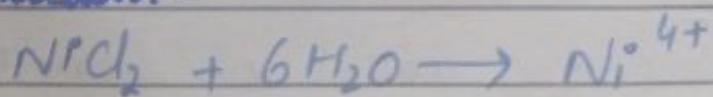

scale + x-axis


$$1 \text{ big box} = 5$$

$$1 \text{ small box} = 0.5$$

y-axis

$$B + B_0 = 0.1$$

$$S \cdot B = 0.01$$

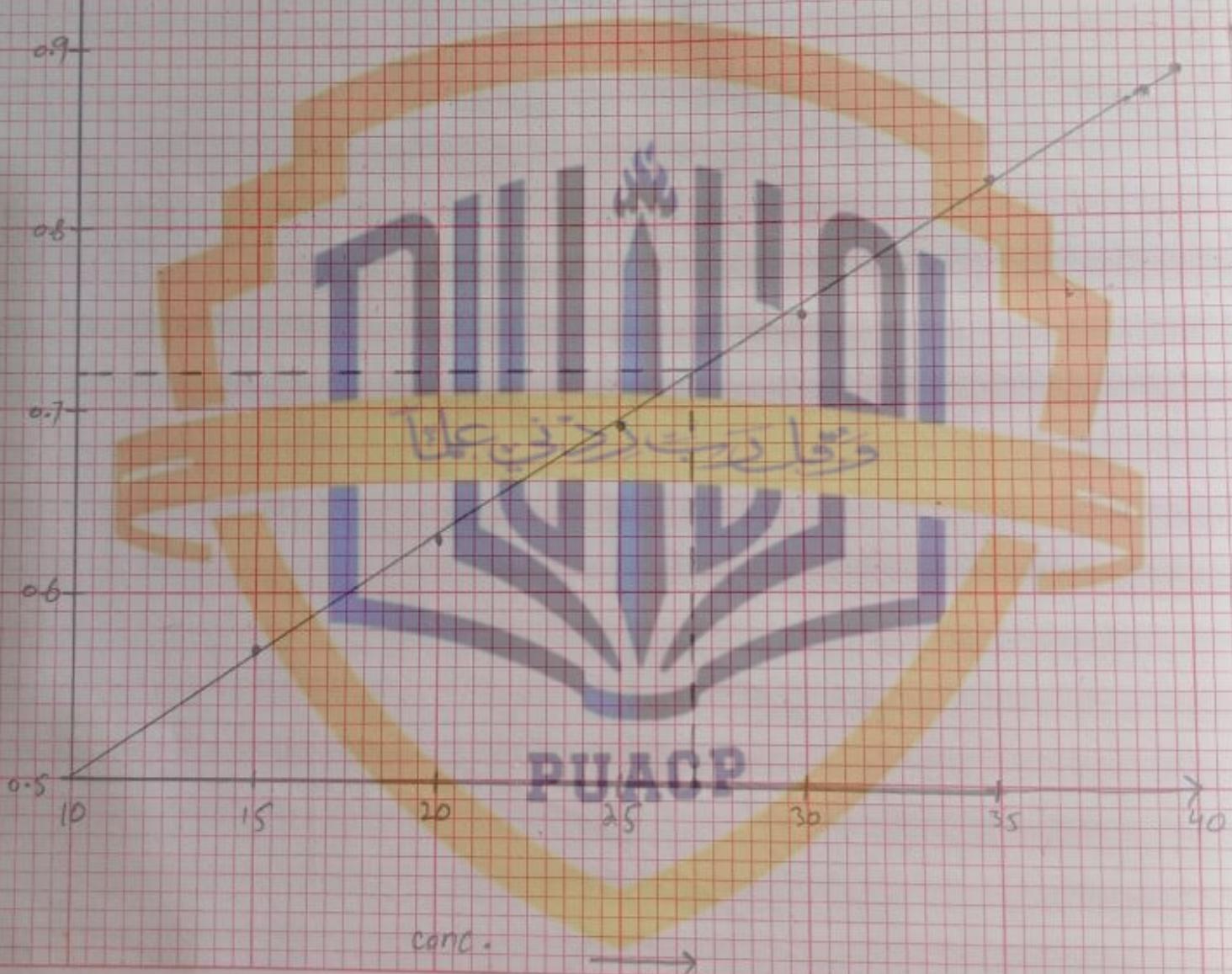


RIZVI NOTE BOOK

## Results:

Conc. of unknown sample is

2ml of liq  $\text{NH}_3$  + 15ml  $\text{H}_2\text{O}$

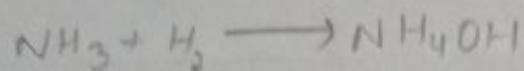


Observation:-

concentration PPm	Absorbance
10	0.5
15	0.56
20	0.629
25	0.695
30	0.755
35	0.820
40	0.885
unknown	0.72

Result:-

conc. of unknown sample is  
27 ppm at absorbance.

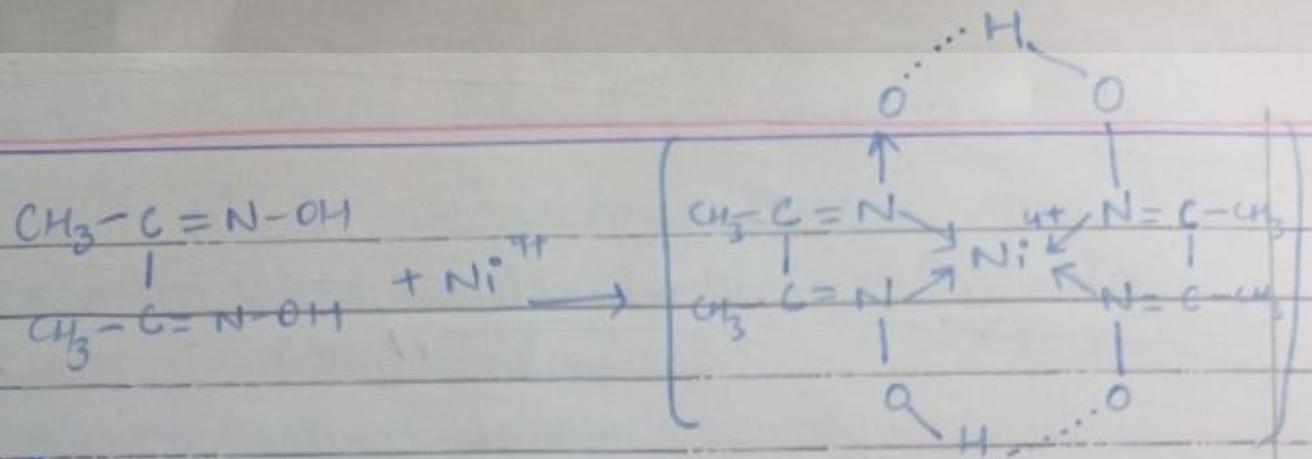

### Procedure:-

- Prepare stock soln of  $\text{NiCl}_2 \cdot 6\text{H}_2\text{O}$ . From this stock solution. Further prepare the conc. soln from 10 ppm - 40 ppm.
- Then taken seven test tubes and tag them.
- Add 5 ml of each ppm conc. of each test tubes, then add 1 ml of ~~Ba~~ water as oxidizing agent. to convert  $\text{Ni}^{2+}$  to  $\text{Ni}^{4+}$
- Then add 1-2 drops of  $\text{NH}_4\text{OH}$  act as a catalyst. Now add 5 ml of complexating agent.
- Allow to stand for 15 min.
- Take readings from spectrophotometer and plot graph.

### Uses of $\text{Ni}^{2+}$

- Chemical element having atomic no 28.
- It is used as metal alloy like steel, brasses and bronze.
- It is used in production of rechargeable batteries, electric generators.

## Experiment # 12

Determine conc. of  $\text{Fe}^{2+}$  by  
colorimetry or spectrophotometry.

### Preparation of reagent

(a) Stock soln:- =  $\frac{\text{Molar mass of compound}}{\text{Molar mass of } \text{Fe}^{2+}}$

$$= \frac{278}{56} = 4.97 \text{ ml H}_2\text{O}$$

(b) Prep of different ppm conc. from stock.

(i) Prep of 10 ppm soln:

Given = Required

$$C_1 V_1 = C_2 V_2$$

$$1000 \times V_1 = 10 \times 100$$

$$V_1 = \frac{1000 \times 10}{100} = 1 \text{ ml / 100 ml H}_2\text{O}$$

(ii) preparation of 15 ml soln:-

$$C_3 V_3 = C_4 V_4$$

$$1000 \times V_1 = 15 \times 100$$

$$V_1 = 1.5 \text{ ml / H}_2\text{O}$$

# Experiment No. 12

Determine concentration of  $\text{Fe}^{2+}$  colorimetrically or spectrophotometrically.

## Theory:-

The measurement of conc. of sample by calibration as a result of standard soln analysis is known as colorimetric analysis.

## Spectroscopy

It is the interaction of electromagnetic radiation with matter.

There are various types of spectroscopy.

## Atomic:-

It is further divided into

- atomic absorption
- atomic emission

## Molecular

Further includes

- UV Visible spectroscopy
- NMR spectroscopy
- IR spectroscopy
- Mass spectroscopy

## Source of excitation

- It can be flamed by hollow cathode lamp.

- Inductively coupled Plasma (ICP)

(3) prep of 20 ppm soln:-

$$C_1 V_1 = C_2 V_2$$

$$1000 \times V_1 = 20 \times 100$$

$$V_1 = 2 \text{ ml} / 100 \text{ ml H}_2\text{O}$$

(4) prep of 25 ppm soln:-

$$C_1 V_1 = C_2 V_2$$

$$1000 \times V_1 = 25 \times 100$$

$$V_1 = 2.5 \text{ ml} / 100 \text{ ml H}_2\text{O}$$

(5) prep of 30 ppm soln:-

$$C_2 V_2 = C_3 V_3$$

$$1000 \times V_2 = 30 \times 100$$

$$V_2 = 3 \text{ ml} / 100 \text{ ml H}_2\text{O}$$

(6) prep of 35 ppm soln

$$C_4 V_4 = C_5 V_5$$

$$0.1 \times 250 = 17 \times V_2$$

$$V_2 = 1.5 \text{ ml}$$

(7) prep of 40 ppm soln

$$C_6 V_6 = C_7 V_7$$

$$1000 \times V_6 = 40 \times 100$$

$$V_6 = 4 \text{ ml} / 100 \text{ ml H}_2\text{O}$$

## Beer Lambert law

It states that


The adsorption of light by sample is directly proportional to conc. and path length of sample.

$$A \propto C \cdot l$$


$$A = a \cdot c \cdot l$$

As we fix path length that is  $l=1\text{ cm}$   
then  $A \propto C$

Calibration Curve

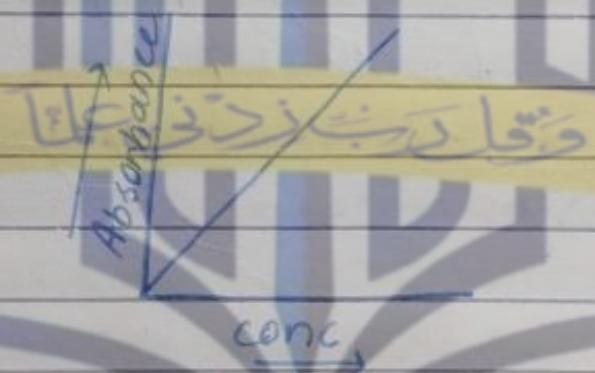

Equation:-

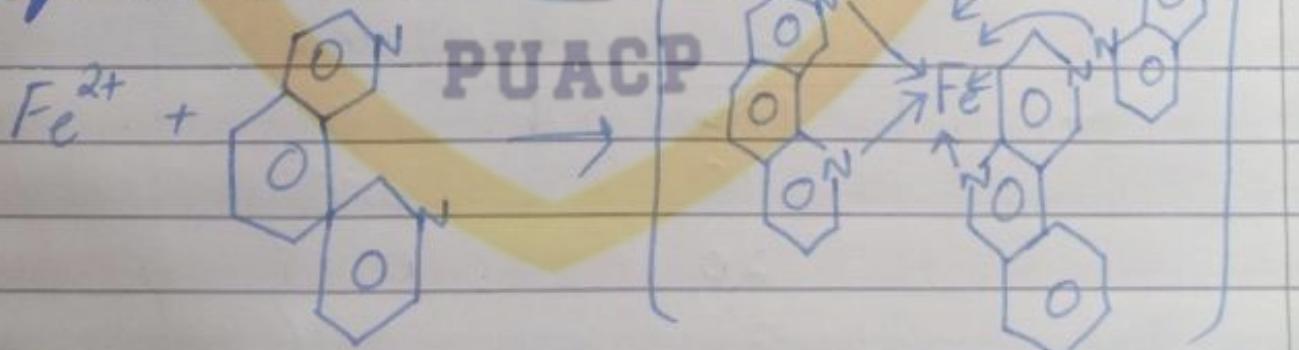

Procedure:-

- Prepare soln of  $\text{FeSO}_4 \cdot 7\text{H}_2\text{O}$  from this Stock soln. Further prepare Seven conc. soln from 10 ppm - 40 ppm.

Scales -  
X-axis

$$1 \text{ B-B} = 5 \text{ ppm}$$

$$1 \text{ small B} = 0.5 \text{ ppm}$$

Y-axis

$$1 \text{ B-B} = 0.3$$

$$1 \text{ s-B} = 0.03$$

Absorbance →

1.5

1.2

0.9

0.6

0.3

0

10

15

20

25

30

35

concentration →

### RIZVI NOTE BOOK

30	1.29
35	1.62
40	1.94
Unknown	0.9

## Preparation of Buffer soln:-

(i) 0.1 M CH<sub>3</sub>COOH :-

$$M_1 V_1 = M_2 V_2$$

$$0.1 \times 250 = 17 \times V_2$$

$$V = 1.8 \text{ ml}$$

(ii) 0.1 M CH<sub>3</sub>COONa

$$0.1 M = n \times 1$$

$$n = 0.82 \text{ g}$$

$\Rightarrow$  Preparation of 1% 8 hydroxyquinoline:-

upto 100 ml with Buffer

$\Rightarrow$  Prep. of PUACP

PPM	Absorbance
10	0.1
15	0.33
20	0.65
25	0.96
30	1.29
35	1.62
40	1.94
Unknown	0.9

- Then take Seven test tubes and tag them.
- Add 5ml of each ppm conc. to each test tubes, then add 1ml of Br-water as oxidizing agent to convert  $\text{Ni}^{2+}$  to  $\text{Ni}^{4+}$ .
- Then add 1-2 drops of  $\text{NH}_4\text{OH}$  act as a catalyst. Now add 5ml of complexating agent i.e 1-10 phenanthroline.
- Allow to stand for 15 min.
- Take reading from spectrophotograph/ spectrophotometer and plot graph.

### Uses of Ferrous

- It is strong metal and used in manufacturing of machine, tools, ships, used in building parts.

### **PUACP**

- Used in manufacturing of tools, surgical equipments and appliances.
- Used in production of steels.

## Experiment # 13:-

Determine conc. of  $\text{Cr}^{3+}$   
colorimetry.

### Preparation of reagent

1. stock soln:  $\frac{\text{molar mass of comp}}{\text{molar mass of}}$

$$= \frac{194.2}{54}$$

$$= 0.373 \text{ g/100 ml}$$

2) Preparation of diff. ppm from stock soln:-

(i) Prep. of 10 ppm soln:-

$$C_1V_1 = C_2V_2$$

$$1000 \times V_1 = 10 \times 100 \quad V_1 = 1 \text{ ml/100 ml}$$

(ii) Prep. of 15 ppm soln:-

$$C_1V_1 = C_2V_2$$

$$1000 \times V_1 = 15 \times 100$$

$$V_1 = 1.5 \text{ ml/100 ml H}_2\text{O}$$

(iii) Prep. of 20 ppm soln:-

$$1000 \times V_1 = 20 \times 100 \text{ ml}$$

$$V_1 = 2 \text{ ml/100 ml H}_2\text{O}$$

# Experiment No. 13

Determine conc. of  $\text{Cr}^{3+}$  colorimetrically.

Theory:-

Measurement of sample by calibration curve are result of standard soln analysis is known as colorimetric analysis. Spectroscopy is the interaction of electromagnetic radiation with matter.

Beer Lambert law


The absorption of light by the sample is directly proportional to conc. and path length.

$$\checkmark A \propto C \cdot l$$

$$A = E \cdot C \cdot l$$

As we know that path length is fixed that  $l = 1\text{ cm}$  then

$$A \propto C$$

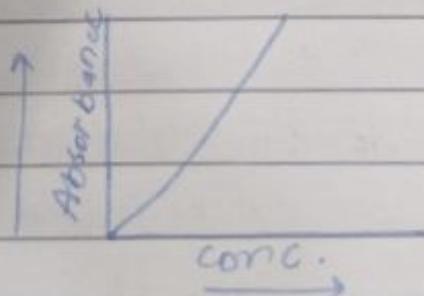

(iv) Prep of 25 ppm:-

$$C_1V_1 = C_2V_2$$

$$V_1 \times 1000 = 25 \times 100$$

$$V_1 = 2.5 \text{ ml}/100 \text{ ml H}_2\text{O}$$

(v) Prep of 30 ppm:-

$$C_1V_1 = C_2V_2$$

$$1000 \times V_1 = 30 \times 100 \text{ ml H}_2\text{O}$$

$$V_1 = 3 \text{ ml}/100 \text{ ml H}_2\text{O}$$

(vi) Prep of 35 ppm:-

$$C_1V_1 = C_2V_2$$

$$1000 \times V_1 = 35 \times 100$$

$$V_1 = 3.5 \text{ ml}/100 \text{ ml H}_2\text{O}$$

(vii) Prep of 40 ppm:-

$$C_1V_1 = C_2V_2$$

$$1000 \times V_1 = 40 \times 100$$

$$V_1 = 4 \text{ ml}/100 \text{ ml H}_2\text{O}$$

Prep. of complexating agent

Take 0.25 ml of diphenyl carbazon  
dissolve it in 30 ml of ethanol and then  
dilute upto 100ml with  $\text{H}_2\text{O}$ .







Draw back:-

spectrophotometer does not give accurate reading of highly concentrated soln.

Source of light:-

In UV/visible spectrophotometer the source is deuterium tungsten lamp.

Equation:-

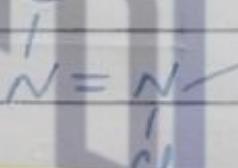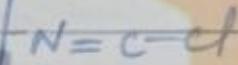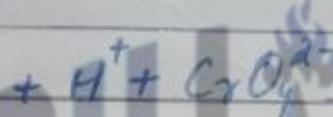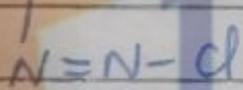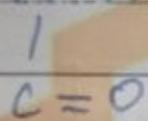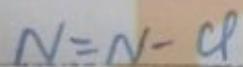

Preparation of soln:-

Given

Required

$$M_1 V_1 = M_2 V_2$$

$$18 V_1 = 6 \times 100$$

$$V_1 = \frac{6 \times 100}{18}$$

$$V_1 = 33.3 \text{ ml / 100 ml H}_2\text{O}$$

Stock solution:-

= Molar mass of compound

Molar mass of K

$$= \frac{194.2}{52} - 3.72 \text{ g / 100 ml H}_2\text{O}$$

scale:- x-axis


$$1 \text{ B.B} = 5 \text{ ppm}$$

$$1 \text{ small B} = 0.5 \text{ ppm}$$

y-axis

$$1 \text{ B.B} = 0.01 \text{ A}$$

$$1 \text{ small B} = 0.001 \text{ A}$$

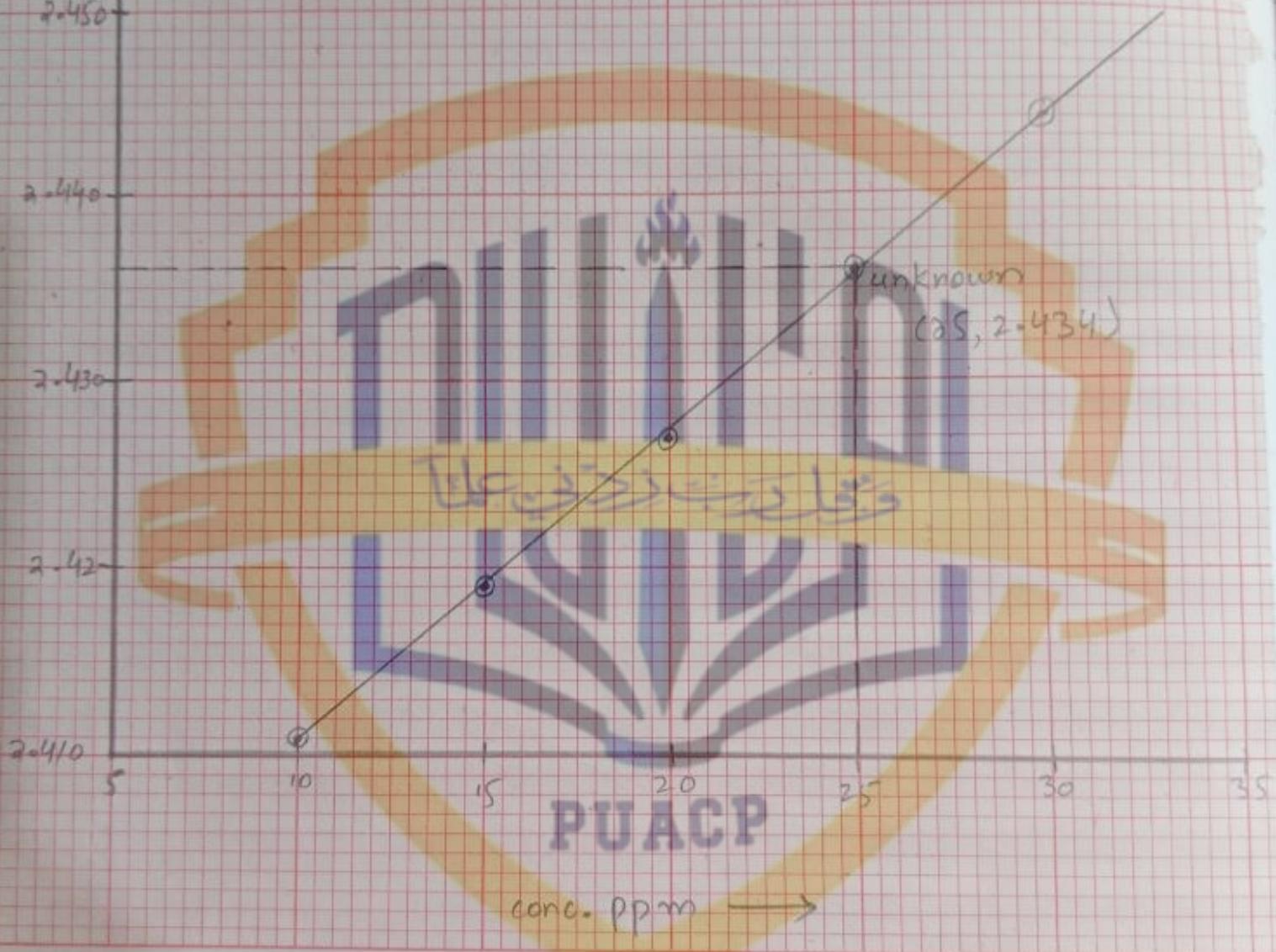

RIZVI NOTE BOOK

Result:-

The conc. of unknown sample  
is 23 ppm at absorbance  
2.430.

prep. of 6M  $H_2SO_4$  :-

$$m_1 v_1 = m_2 v_2$$

$$18 \times v_1 = 6 \times 100$$

$$v_1 = \frac{6 \times 100}{18}$$

$$v_1 = 33.3 \text{ ml / 100 ml } H_2O.$$

Observations:-

Conc.	Absorbance
10 جنی علی	2.470
15	2.418
20	2.425
25	2.433
30	2.441
35	2.499
40	2.458
unknown.	2.430

Result:-

The conc. of unknown sample  
is 23 ppm at absorbance  
2.430.

$$= 0.373/100 \text{ ml} = 1000 \text{ PPm.}$$

From this stock soln further prepare seven conc. from 10ppm-40ppm.

## Preparation of complexating agent

Take 0.25 g of diphenyl carbazole  
dissolve it into 30ml of ethanol and  
then dilute upto 100ml of H<sub>2</sub>O.

### Procedure:-

- Prepare a stock soln of K<sub>2</sub>C<sub>2</sub>O<sub>4</sub> from this stock soln prepare seven conc. soln from 10 ppm - 40 ppm.
- Then take 7 test tubes and add 5 ml of each ppm soln in each tube. Now add 1 ml of 6 M H<sub>2</sub>SO<sub>4</sub> and 5 ml of complexing agent i.e. diphenyl carbazole.
- Allow it to stand for 15 min.
- Take reading or check absorbance in spectrophotometer or calorimeter and then plot a calibration curve.  
→ The color of complex is violet.

## ⇒ Uses of Cr

- Used to harden the steel in manufacturing of stainless steel.
- used to produce alloy
- Cr-plating used to give a polish mirror finish to steel.
- Cr-plate ~~that is used in~~ bathroom fitting.
- Cr-plated car and lorries parts are also common such as bumper.

## Experiment # 14:-

Determine amount / l of  $\text{Cu}^{2+}$  gravimetrically.

1) preparation of Reagent:-

$\Rightarrow 1\% \text{ CuSO}_4$  soln:-

Dissolve 1g of  $\text{CuSO}_4$  in  $\text{H}_2\text{O}$  and make volume upto 100ml.

$\Rightarrow 1\% \text{ Salicyldoxime}$

Dissolve 1g of salicyldoxime in 10.15 ml ethanol for completely soluble and make vol. upto 100ml.

Calculations:-

$$\text{wt. of filter paper} = w_1 = 1.21$$

$$\text{wt. of filter paper} = w_2 = 1.27$$

$$\begin{aligned}\text{P.U.A.C.P.} &= w_2 - w_1 \\ &= 1.27 - 1.21 \\ &= 0.15 \text{ g}\end{aligned}$$

335g of complex contain  $\text{Cu}^{2+} = 64$

$$1g \quad " \quad " \quad " \quad = \frac{64}{335}$$

## Experiment No. 14

$\text{Cu}^{2+}$  Determine amount /l of gravimetrically.

### Theory:-

In this technique amount of soluble species is determined by converting it into insoluble species by using some precipitating agents which may be organic or inorganic reagents.

The amount of soluble species is calculated by using the knowledge of molar weight of species formed.

### Use of gravimetric technique

It is used on wide scale in fertilizer industry ~~for determination of nutrients~~ PUAC

In cement industry for determination of silica,  $\text{Al}^{3+}$ ,  $\text{Ca}^{2+}$

In electroplating for determination of  $\text{Zn}^{2+}$ ,  $\text{Ca}^{2+}$ ,  $\text{Ag}^+$  etc.

$$= 0.15$$

$$0.15 \text{ g} \quad " \quad = \frac{64}{335} \times 0.15$$

20 ml of sample contain  $\text{Cu}^{2+} = 0.028 \text{ g}$


Result:-

$$\begin{aligned} 1 \text{ ml} & \quad " \quad = \frac{0.028}{20} \times 1000 \\ & = 1.49 \text{ g/l} \end{aligned}$$

The amount /l of  $\text{Cu}^{2+}$  in  
given Sample is 1.49 g/l.

PUACP

Equation:-

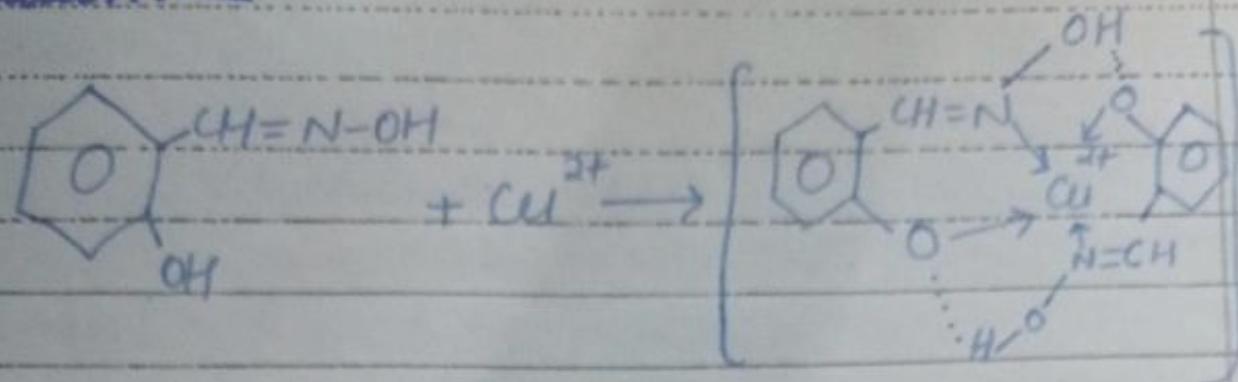

Preparation of soln:-

1%  $\text{CuSO}_4$  soln

Dissolve 1 g of  $\text{CuSO}_4$  in  $\text{H}_2\text{O}$  and make the volume upto 100 ml.

1% Salicylaldoxime

Dissolve 1 g of salicylaldoxime in 10-15 ml ethanol per the complete solubility of it and make volume upto 100ml of  $\text{H}_2\text{O}$ .

Procedure:- PUACP

- Take 20 ml of  $\text{CuSO}_4$  soln in tank.
- Add few drops of Salicylaldoxime soln until ppt formed.
- Heat soln so that precipitates are coagulated. Filter it and yield is weight of ppt.

- ppt are light green yield of ppt is 0.15g

### Uses of $\text{Cu}^{2+}$

- In building constructions used for power generation and transmission.
- Used for manufacturing of electric power.
- used in production of industrial machinery and transportation vehicles.