

Chapter 23 - Potentiometry

- Read: pp. 659-679 Problems: 23-2,4,7,13,14
- Potentiometric methods are based upon measurements of the **potential** of electrochemical cells in the **absence** of appreciable currents (an equilibrium measurement, therefore, the Nernst equation is applicable).
- All equipment is simple: an **indicator** electrode, a **reference** electrode and a **potential measuring device**.
- Billions of these measurements are made annually. Importance in environmental and medical applications. For example, pH, ion selective electrodes, blood gas analysis (O_2 , CO_2), etc.

Reference Electrodes

Ideal Properties

1. Reversible and obeys Nerst eq.
2. Stable potential with time.
3. Returns to original position after passage of small currents.
4. Little hysteresis with temperature.

$\text{Hg}/\text{Hg}_2\text{Cl}_2(\text{sat'd}), \text{KCl} (x\text{M}) //$
 $E^\circ = 0.244 \text{ V vs. NHE}$

$\text{Ag}/\text{AgCl}(\text{sat'd}), \text{KCl} (x\text{M}) //$
 $E^\circ = 0.222 \text{ V vs. NHE}$

Know your redox reactions!!

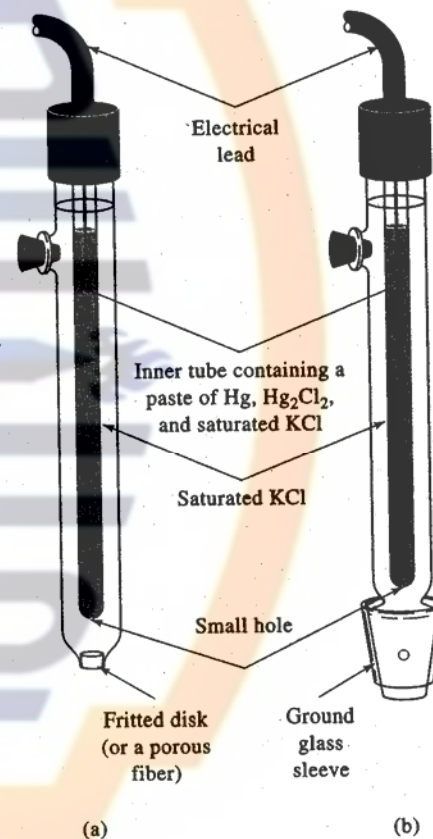
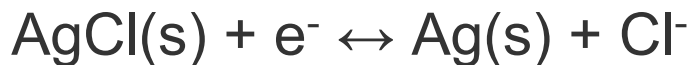


Figure 23-1 Typical commercial calomel reference electrodes.


Reference is considered the anode!

Types of Metallic Indicator Electrodes

Metallic and Membrane


Electrodes of a First Kind

metallic electrodes in direct equilibrium with the cation derived from the metal.

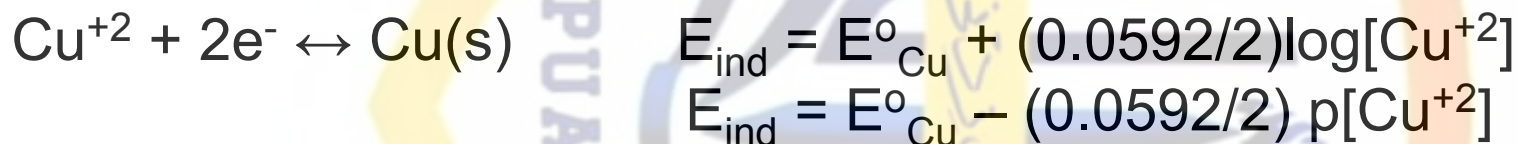

Electrodes of a Second Kind

metallic electrode that is responsive to the activity of an anion with which its ion forms a precipitate or stable complex ion.

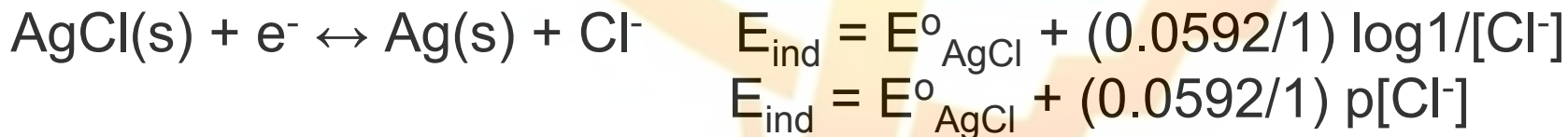

Remember that inverting the log term changes the sign in front of it.

$$\text{pX} = -\log X$$

Types of Metallic Indicator Electrodes

Electrodes of a Third Kind

A metal electrode, under certain circumstances can be made to respond to a different cation.



$$K_f = (\text{aCa}^{+2} \cdot \text{aY}^{-4})/\text{aCaY}^{-2}$$

Metallic Indicator Electrodes

Pt, Au, Pd (inert metals) are responsive to the activities of the oxidized and reduced forms of the redox couple near the electrode surface.




$$E = E^{\circ} + (0.0592/n) \log ([\text{Ce}^{+4}]/[\text{Ce}^{+3}])$$

The pH Electrode – Membrane Electrode

TABLE 23-2 Types of Ion-Selective Membrane Electrodes

- | |
|--|
| A. Crystalline Membrane Electrodes |
| 1. Single crystal
Example: LaF_3 for F^- |
| 2. Polycrystalline or mixed crystal
Example: Ag_2S for S^{2-} and Ag^+ |
| B. Noncrystalline Membrane Electrodes |
| 1. Glass
Examples: silicate glasses for Na^+ and H^+ |
| 2. Liquid
Examples: liquid ion exchangers for Ca^{2+} and neutral carriers for K^+ |
| 3. Immobilized liquid in a rigid polymer
Examples: polyvinyl chloride matrix for Ca^{2+} and NO_3^- |

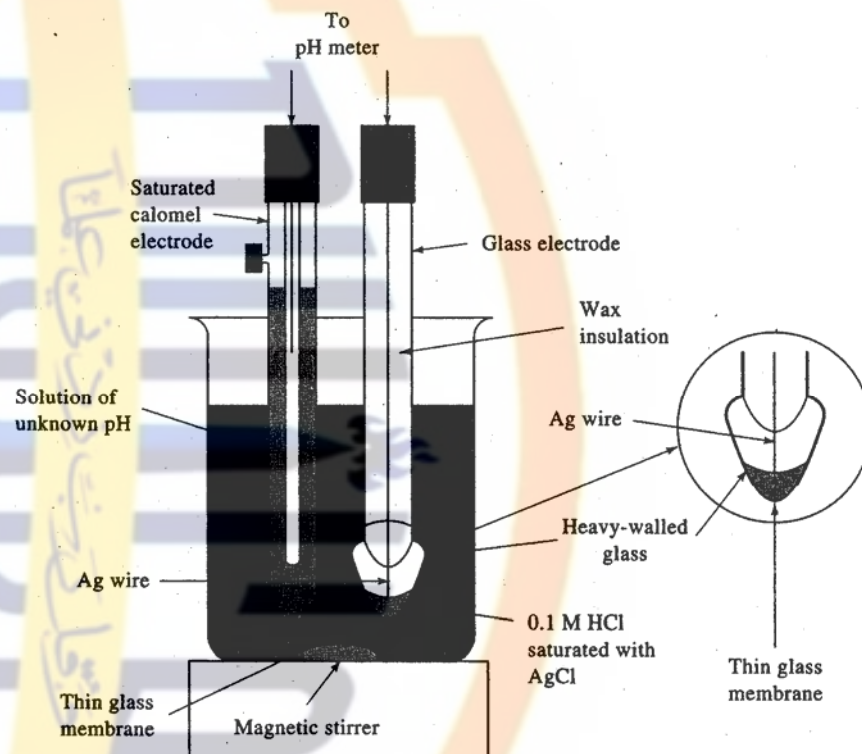

The membranes must have (i) minimal solubility, (ii) be ion conductors and not electrical conductors, and (iii) some selective interaction with the analyte of interest.

Membranes Selectively Separate Charge

Charge separation = potential!

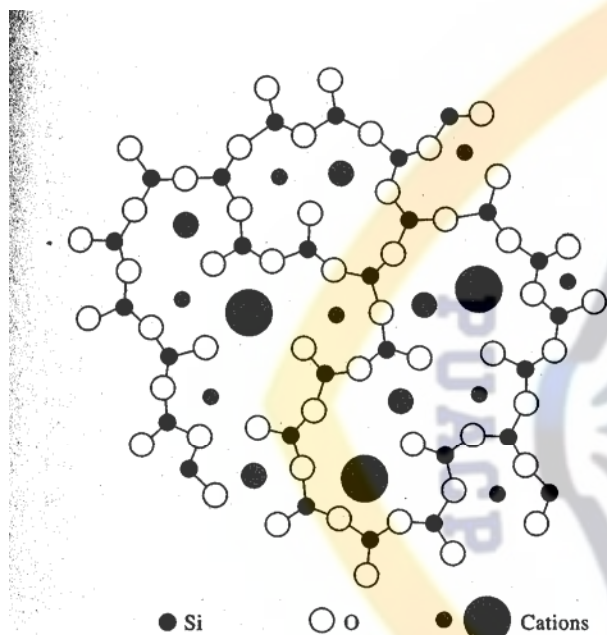


Figure 23-5 Cross-sectional view of a silicate glass structure. In addition to the three Si—O bonds shown, each silicon is bonded to an additional oxygen atom, either above or below the plane of the paper. (Adapted with permission from G. A. Perley, Anal. Chem., 1949, 21, 395. Copyright 1949 American Chemical Society.)

Hygroscopic gel layers



$$E_b = E_1 - E_2 = 0.0592 \log a_1/a_2$$

$$E_b = L' + 0.0592 \log a_1$$

$$= L' - 0.0592 \text{ pH}$$

$$\text{where } L' = -0.0592 \log a_2$$




The pH Electrode

- Alkaline error – sensitive to alkali metal ions at pH greater than 12.
- Acid error – at pH less than 0.5, values obtained with the pH electrode are high.
- Dehydration – must keep membrane moist.
- Errors in low ionic strength – varying junction potentials.
- Errors in pH of standard buffer solutions.

Reference Electrode 1 // $[H_3O^+] = a_1$ / membrane / $[H_3O^+] = a_2$ // Reference Electrode 2

Boundary potential, E_b , is sensitive to solution pH!

Crystalline Membrane Electrodes

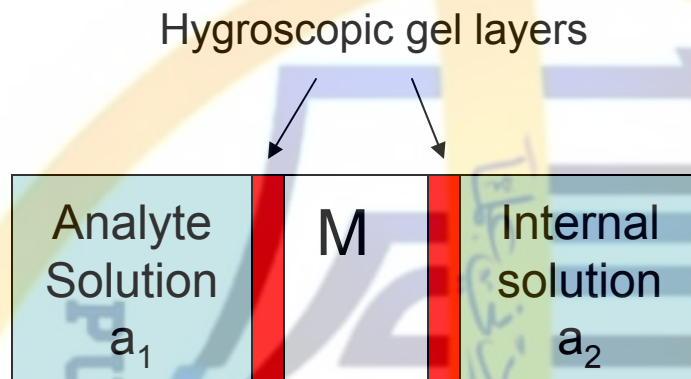

$$E_b = E_1 - E_2 = 0.0592 \log a_1/a_2$$

$$E_{\text{ind}} = L - 0.0592 \log a_{\text{F}^-} = L + 0.0592 \text{p[F]}$$

Crystalline Membrane Electrodes

TABLE 23-3 Commercial Solid-State Electrodes^a

Analyte Ion	Concentration Range, M	Interferences ^b
Br ⁻	10 ⁰ to 5 × 10 ⁻⁶	mr: 8 × 10 ⁻⁵ CN ⁻ ; 2 × 10 ⁻⁴ I ⁻ ; 2 NH ₃ ; 400 Cl ⁻ ; 3 × 10 ⁴ OH ⁻ . mba: S ²⁻
Cd ²⁺	10 ⁻¹ to 10 ⁻⁷	Fe ²⁺ + Pb ²⁺ may interfere. mba: Hg ²⁺ , Ag ⁺ , Cu ²⁺
Cl ⁻	10 ⁰ to 5 × 10 ⁻⁵	mr: 2 × 10 ⁻⁷ CN ⁻ ; 5 × 10 ⁻⁷ I ⁻ ; 3 × 10 ⁻³ Br ⁻ ; 10 ⁻² S ₂ O ₃ ²⁻ ; 0.12 NH ₃ ; 80 OH ⁻ . mba: S ²⁻
Cu ²⁺	10 ⁻¹ to 10 ⁻⁸	high levels Fe ²⁺ , Cd ²⁺ , Br ⁻ , Cl ⁻ . mba: Hg ²⁺ , Ag ⁺ , Cu ⁺
CN ⁻	10 ⁻² to 10 ⁻⁶	mr: 10 ⁻¹ I ⁻ ; 5 × 10 ³ Br ⁻ ; 10 ⁶ Cl ⁻ . mba: S ²⁻
F ⁻	sat'd to 10 ⁻⁶	0.1 M OH ⁻ gives <10% interference when [F ⁻] = 10 ⁻³ M
I ⁻	10 ⁰ to 5 × 10 ⁻⁸	mr: 0.4 CN ⁻ ; 5 × 10 ³ Br ⁻ ; 10 ⁵ S ₂ O ₃ ²⁻ ; 10 ⁶ Cl ⁻
Pb ²⁺	10 ⁻¹ to 10 ⁻⁶	mba: Hg ²⁺ , Ag ⁺ , Cu ²⁺
Ag ⁺ /S ²⁻	10 ⁰ to 10 ⁻⁷ Ag ⁺ 10 ⁰ to 10 ⁻⁷ S ²⁻	Hg ²⁺ must be less than 10 ⁻⁷ M
SCN ⁻	10 ⁰ to 5 × 10 ⁻⁶	mr: 10 ⁻⁶ I ⁻ ; 3 × 10 ⁻³ Br ⁻ ; 7 × 10 ⁻³ CN ⁻ ; 0.13 S ₂ O ₃ ²⁻ ; 20 Cl ⁻ ; 100 OH ⁻ . mba: S ²⁻

^aFrom: *Handbook of Electrode Technology*, pp. 10-13, Appendix, Orion Research: Cambridge, MA, 1982. With permission.

^bmr: maximum ratio $\left(\frac{C_{\text{interferent}}}{C_{\text{analyte}}}\right)$ for no interference.

mba: must be absent.

Liquid Membrane Electrodes


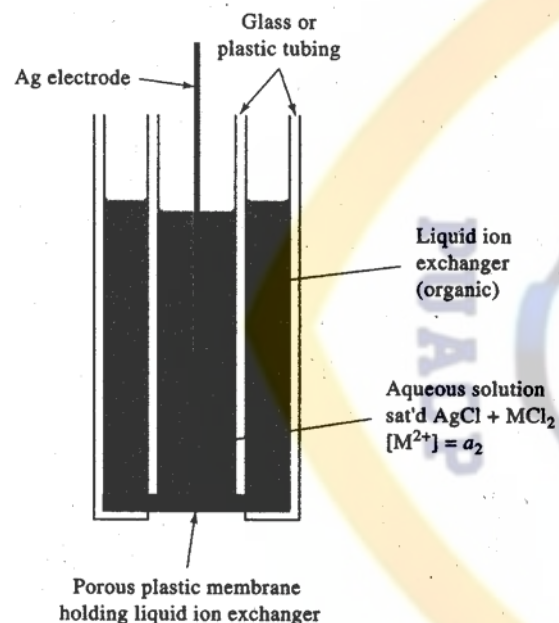


Figure 23-8 Liquid membrane electrode sensitive to M²⁺.

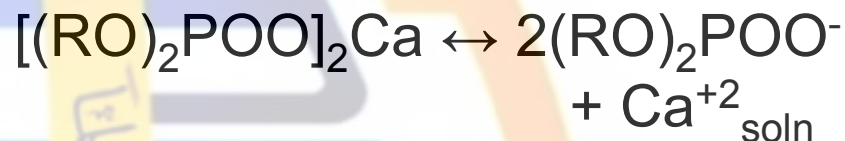

Fill membrane with a compound that selectively binds with the analyte.

$$E_{ind} = L + (0.0592/2) \log a_1$$

(Ca⁺² is divalent)

$$E_{ind} = L' - (0.0592/2) p[Ca]$$



Gas Sensing Electrodes

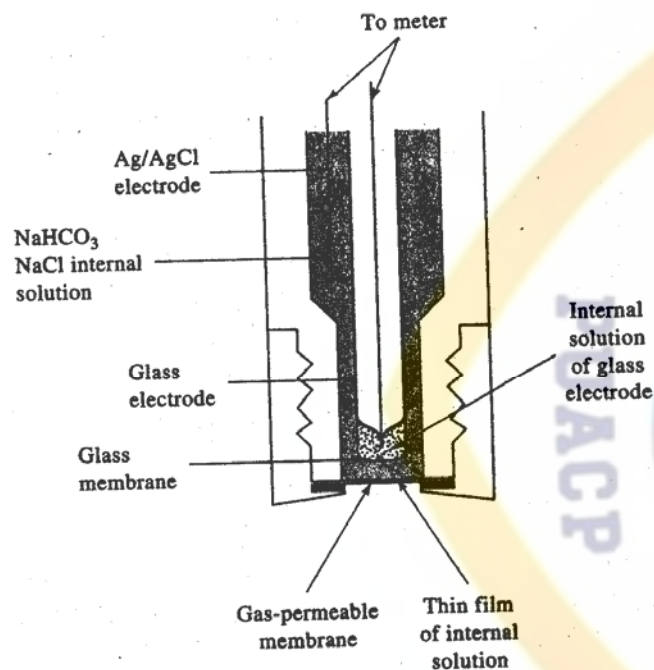


TABLE 23-5 Commercial Gas-Sensing Probes

Gas	Equilibrium in Internal Solution	Sensing Electrode
NH ₃	$\text{NH}_3 + \text{H}_2\text{O} \rightleftharpoons \text{NH}_4^+ + \text{OH}^-$	Glass, pH
CO ₂	$\text{CO}_2 + \text{H}_2\text{O} \rightleftharpoons \text{HCO}_3^- + \text{H}^+$	Glass, pH
HCN	$\text{HCN} \rightleftharpoons \text{H}^+ + \text{CN}^-$	Ag ₂ S, pCN
HF	$\text{HF} \rightleftharpoons \text{H}^+ + \text{F}^-$	LaF ₃ , pH
H ₂ S	$\text{H}_2\text{S} \rightleftharpoons 2\text{H}^+ + \text{S}^{2-}$	Ag ₂ S, pS
SO ₂	$\text{SO}_2 + \text{H}_2\text{O} \rightleftharpoons \text{HSO}_3^- + \text{H}^+$	Glass, pH
NO ₂	$2\text{NO}_2 + \text{H}_2\text{O} \rightleftharpoons \text{NO}_2^- + \text{NO}_3^- + 2\text{H}^+$	Immobilized ion exchange, pNO ₃

Figure 23-11 Schematic of a gas-sensing probe for carbon dioxide.